



Emergency Medical Systems





INDEX

	PAGE
MEDICAL SUCTION UNIT	6
EMERGENCY BAGS & RUCKSACKS	24
IMMOBILIZATION	39
RESUSCITATOR BAGS	56
DIAGNOSTICS	63
REANIMATION	73
MEDICATION	83
MEDICAL OXYGEN	127
OTHERS, ACCESSORIES, DIDACTICS	134
CONTACTS AND SELLING CONDITIONS	139

Dear Customer,
herewith we inform you that due to our Ministry of Health instructions related to the advertisement of medical devices, diagnostics in vitro and surgical medical devices, the featuring of the pages hereafter are appointed to professionals. We also remind you that our products description only aims to inform and are not intended for promotional purpose.

INTRODUCTION



Oscar Boscarol srl has always encouraged the research for vanguard solutions and technologies with compelling quality that might match the operators needs during emergency activities. This is the force that led our products to become a helping hand, an indispensable support and constant reference for all insiders. Today the Oscar Boscarol srl is a worldwide leader exporting half of its products to foreign markets.





MADE IN ITALY

Our OB medical suction units are entirely manufactured in Italy at our headquarters in Bolzano in order to guarantee care and quality control, starting from the materials used up to the assembly of the device and the control of its correct functioning before customer delivery.

HISTORY-RELIABILITY

We have been manufacturing medical suction units since the 90s and we are at the top of the international competition with over 40,000 units sold worldwide (Europe, Asia, India, South America and Africa).

COMPLIANCE WITH REGULATIONS

All our devices comply with current regulations: our quality system was developed according to the ISO 9001 and ISO 13485 standards, by always placing first the need to produce qualitatively impeccable devices that comply with the mandatory requirements.

DIRECTIVES, RULES AND REGULATIONS OF REFERENCE

Current regulatory references our suction units comply with.

MDD 93/42/CEE	Council Directive concerning medical devices
ISO 10079-1	Medical suction equipment - Part 1: electrically powered suction equipment
IEC 60601-1	Medical Electrical Equipment. General requirements for basic safety and essential performance
IEC 60601-1-2	Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral Standard: Electromagnetic disturbances - Requirements and tests
DO-160	Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral Standard: Electromagnetic disturbances - Requirements and tests ONLY FOR OUR AVIO SUCTION UNIT MODELS
IEC 60601-1-11	Medical Electrical Equipment General requirements for basic safety and essential performance - Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment
IEC 60601-1-12	Medical Electrical Equipment General requirements for basic safety and essential performance - Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the emergency medical services environment
EN 1789	Medical vehicles and their equipment - Road ambulances
ECE R10 REV. 5	Uniform provisions concerning the approval of vehicles with regard to electromagnetic compatibility



OB 3000

OB 3000

FLOW RATE-MAX VACUUM	33 LPM - 800 mbar (80 kPa)
RATED VOLTAGE	11÷30 Vcc
AUTONOMY	60 minutes ±10%
IP RATING	IP 44
BATTERY TYPE	LIPO 11,1V - 5 Ah removable
WARRANTY	5 years
PERIODIC MAINTENANCE	2 years (automatic reminder)
RTCA D0160G	YES in AVIO version
DIMENSIONS	360×244(h)×110 mm
WEIGHT	2,6 kg with accessories (without wall bracket)



OB 3000

POWERFUL AND RELIABLE

Max. flow rate 33 liter/minute with double piston suction pump and max. vacuum rate 800 mbar. Inner microprocessor records production data and functioning tracing and a safety inspection reminder indicates that the suction unit has to be sent for the scheduled maintenance. Analogic fluorescent vacuum meter, with mmHg scale upon request. Analogic fluorescent vacuum meter, with mmHg scale upon request.



BATTERY AND RECHARGING SYSTEM

Battery can be easily replaced by user and battery status indicator with no need to turn the device on.

Upon request the device can be supplied with a wall bracket (conformed to EN 1789) and with a 2-poles main power supply male connector.

WEIGHT AND DIMENSIONS

FA version: 1000 ml autoclavable canister with protection overflow valve and protection filter directly on the cover.

FM version: 1000 ml OB-J canister suitable for the use of Serres liners. Extreme lightwightness and reduced dimensions make of this unit an easily portable device.



BSU3000	- OB 3000 FA MEDICAL SUCTION UNIT
BSU3020	- OB 3000 FM MEDICAL SUCTION UNIT
BSU3000A	- OB 3000 AVIO FA MEDICAL SUCTION UNIT
BSU3020A	- OB 3000 AVIO FM MEDICAL SUCTION UNIT



OB MINI

FLOW RATE-MAX VACUUM	25 LPM - 850 mbar (85 kPa)
RATED VOLTAGE	11÷30 Vcc
AUTONOMY	60 minutes ±10%
IP RATING	IP 43
BATTERY TYPE	LIPO 11,1V - 5 Ah
WARRANTY	3 years
PERIODIC MAINTENANCE	2 years (automatic reminder)
RTCA D0160G	YES in AVIO version
DIMENSIONS	290×240(h)×95 mm
WEIGHT	1,8 kg with accessories



OB MINI

RELIABLE AND ULTRA LIGHTWEIGHT

Max. flow rate 25 liter/minute with double piston suction pump and max. vacuum rate 850 mbar. Inner microprocessor records production data and functioning tracing and a safety inspection reminder indicates that the suction unit has to be sent for the scheduled maintenance. Only 1,8 kg for a device complete with 500cc canister: maximum lightweight and compactness.



JAR

500 cc autoclavable canister or upon request 1 l FA and 1 l FM version.
FA version: 1000 ml autoclavable canister with protection overflow valve and protection filter directly on the cover.
FM version: 1000 ml OB-J canister



POWER SUPPLY AND BAG

2-poles main power supply
male connector.
Blue transport bag with
lateral poche with velcro.

BSU360EU - OB MINI 500ML MEDICAL SUCTION UNIT EU
BSU360JP - OB MINI 500ML MEDICAL SUCTION UNIT JP/USA
BSU360UK - OB MINI 500ML MEDICAL SUCTION UNIT UK
BSU370 - OB MINI AVIO 500ML MEDICAL SUCTION UNIT EU
BSU370JP - OB MINI AVIO 500ML MEDICAL SUCTION UNIT JP
BSU370UK - OB MINI AVIO 500ML MEDICAL SUCTION UNIT UK
BSU364EU - OB MINI FA MEDICAL SUCTION UNIT EU
BSU364JP - OB MINI FA MEDICAL SUCTION UNIT JP
BSU364UK - OB MINI FA MEDICAL SUCTION UNIT UK
BSU368EU - OB MINI FM MEDICAL SUCTION UNIT EU
BSU368JP - OB MINI FM MEDICAL SUCTION UNIT JP
BSU368UK - OB MINI FM MEDICAL SUCTION UNIT UK
BSU372 - OB MINI AVIO FA MEDICAL SUCTION UNIT EU
BSU372JP - OB MINI AVIO FA MEDICAL SUCTION UNIT JP
BSU372UK - OB MINI AVIO FA MEDICAL SUCTION UNIT UK
BSU374 - OB MINI AVIO FM MEDICAL SUCTION UNIT EU
BSU374JP - OB MINI AVIO FM MEDICAL SUCTION UNIT JP
BSU374UK - OB MINI AVIO FM MEDICAL SUCTION UNIT UK



OB 2012

OB2012

FLOW RATE-MAX VACUUM	30 LPM - 800 mbar (80 kPa)
RATED VOLTAGE	12÷15 Vcc
AUTONOMY	60 minutes ±10%
IP RATING	IP 34d
BATTERY TYPE	SLA 12V - 4 Ah
WARRANTY	3 years
PERIODIC MAINTENANCE	2 years
RTCA DO 160G	NO
DIMENSIONS	350×240(h)×120 mm
WEIGHT	4,6 kg with accessories (without wall bracket)



HARD AND RESISTANT

Max. flow rate 30 liter/minute with double piston suction pump and max. vacuum rate 800 mbar.



FLACONI

FA version: 1000 ml autoclavable canister with protection overflow valve and protection filter directly on the cover.
FM version: 1000 ml OB-J canister suitable for the use of Serres liners.



RECHARGING SYSTEM

Upon request the device can be supplied with a wall bracket (conformed to EN 1789) and with a 2-poles main power supply male connector.
Red transport bag.

BSU150	- OB2012 FA MEDICAL SUCTION UNIT
BSU154	- OB2012 FA MEDICAL SUCTION UNIT + WALL BRACKET
BSU158	- OB2012 FA MEDICAL SUCTION UNIT COMPLETE
BSU100	- OB2012 LINER MEDICAL SUCTION UNIT
BSU104	- OB2012 LINER MEDICAL SUCTION UNIT + WALL BRACKET
BSU108	- OB2012 LINER MEDICAL SUCTION UNIT COMPLETE

Codes content summary table on p. 17



OB 1000

OB1000

FLOW RATE-MAX VACUUM	25 LPM - 850 mbar (85 kPa)
RATED VOLTAGE	12÷15 Vcc
AUTONOMY	60 minutes ±10%
AUTONOMY	IP 34d
BATTERY TYPE	SLA 12V - 4 Ah
WARRANTY	2 years
PERIODIC MAINTENANCE	2 years
RTCA D0160G	YES in AVIO version
DIMENSIONS	320×250(h)×100 mm
WEIGHT	3,5 kg with accessories (without wall bracket)



RELIABLE AND AVIO VERSION

Max. flow rate 25 liter/minute with double piston suction pump and max. vacuum rate 850 mbar. Available in AVIO version: conformed to RTCA DO 160G.



JAR

FA version: 1000 ml autoclavable canister with protection overflow valve and protection filter directly on the cover.
FM version: 1000 ml OB-J canister suitable for the use of Serres liners.



RECHARGING SYSTEM

Upon request the device can be supplied with a wall bracket (conformed to EN 1789) and with a 2-poles main power supply male connector.
Red transport bag.

BSU210	- OB1000 FA MEDICAL SUCTION UNIT
BSU216	- OB1000 FA MEDICAL SUCTION UNIT EU
BSU220	- OB1000 FA MEDICAL SUCTION UNIT COMPLETE
BSU232	- OB1000 AVIO FA MEDICAL SUCTION UNIT
BSU212	- OB1000 FM MEDICAL SUCTION UNIT
BSU226	- OB1000 FM MEDICAL SUCTION UNIT EU
BSU228	- OB1000 FM MEDICAL SUCTION UNIT COMPLETE
BSU230	- OB1000 AVIO FM MEDICAL SUCTION UNIT

Codes content summary table on p. 16

OB1000 CODES INDEX



OB1000



12 VCC CABLE

BSU855: 2 poles
recharging supplying
cable



WALL BRACKET

BSU810: wall bracket



ALIMENTATORE

BSU895EU: 2 poles
power supply 14,0 V

BSU210

OB1000 FA MEDICAL
SUCTION UNIT



-

-

BSU216

OB1000 FA MEDICAL
SUCTION UNIT EU

-



UK and JP/USA connections
available upon request

BSU220

OB1000 FA MEDICAL
SUCTION UNIT COMPLETE



-

BSU232

OB1000 AVIO FA MEDICAL
SUCTION UNIT



-

-

BSU212

OB1000 FM MEDICAL
SUCTION UNIT



-

-

BSU226

OB1000 FM MEDICAL
SUCTION UNIT EU

-



UK and JP/USA connections
available upon request

BSU228

OB1000 FM MEDICAL
SUCTION UNIT COMPLETE



-

BSU230

OB1000 AVIO FM MEDICAL
SUCTION UNIT



-

-

IT IS POSSIBLE TO REQUEST MAIN POWER SUPPLY ADAPTERS WITH PLUGS FOR UK AND JP/USA SOCKETS

OB2012 CODES INDEX



OB2012



12 VCC CABLE

BSU855: 2 poles
recharging supplying
cable



WALL BRACKET

BSU810: wall bracket



POWER SUPPLY

BSU895EU: 2 poles
power supply 14,0 V

BSU150

OB2012 FA MEDICAL
SUCTION UNIT



-

-

BSU154

OB2012 FA MEDICAL
SUCTION UNIT + WALL
BRACKET



-

BSU158

OB2012 FA MEDICAL
SUCTION UNIT COMPLETE



UK and JP/USA connections
available upon request

BSU100

OB2012 LINER MEDICAL
SUCTION UNIT



-

-

BSU104

OB2012 LINER MEDICAL
SUCTION UNIT + WALL
BRACKET



-

BSU108

OB2012 LINER MEDICAL
SUCTION UNIT COMPLETE



UK and JP/USA connections
available upon request

IT IS POSSIBLE TO REQUEST MAIN POWER SUPPLY ADAPTERS WITH PLUGS FOR UK AND JP/USA SOCKETS

OB500 - STATIONARY MEDICAL SUCTION UNIT

JAR

FA version: 1000 ml autoclavable canister with protection overflow valve and protection filter directly on the cover.
FM version: 1000 ml OB-J canister suitable for the use of Serres liners.

STATIONARY SUCTION UNIT FOR AMBULANCE

Max. flow rate 30 liter/minute and max. vacuum rate 800 mbar.

Dimensions: of suction unit: 175x175x100 mm

Weight: max. 2 kg including mounting parts.



MODULE OR EMBEDDING

Vacuum control and regulation unit available in EMBEDDED version to be embedded on the wall or in MODULE version with ABS module.

-
- BSU464 - OB500 FM MEDICAL SUCTION UNIT UNDERWALL EMBEDDING
 - BSU442 - OB500 FA MEDICAL SUCTION UNIT UNDERWALL EMBEDDING
 - BSU462 - OB500 LINER MEDICAL SUCTION UNIT UNDERWALL EMBEDDING
 - BSU414 - OB500 FM MEDICAL SUCTION UNIT MODULE
 - BSU402 - OB500 FA MEDICAL SUCTION UNIT MODULE
 - BSU412 - OB500 LINER MEDICAL SUCTION UNIT MODULE
-

JARS FOR BOSCAROL SUCTION UNIT



OB-J CANISTER FOR SERRES LINERS - BSU506

New canister OB-J for secretion collection made in polypropylene for medical use. The jar features one canister and one adapter that allows the use of the jar with the disposable liner produced by the company Serres (non included). Complete with "L" joint. The canister is vapour autoclavable at a max temperature of 121°C and at a pressure of 2bar (200kPa). The jar can also be disinfected accordingly to local hospital procedures. Capacity: 1000 ml.



OB-J FA CANISTER - BSU500

New autoclavable suction jar (121° C - 2 bar max.) with made in polypropylene for medical use. Complete of overflow security valve. The protection filter and the angular joint connector for the aspiration tube patient side are directly connected to the cover. Capacity of 1000 ml.



500 ML CANISTER - BSU504

Suction jar for medical suction units, made in autoclavable and shockproof Macrolon with pressure closing lid made in silicone-rubber. Complete with overflow valve. Capacity: 500 ml.

ACCESSORIES FOR BOSCAROL SUCTION UNIT



ASPIRATION PROBE

Jankauer rigid aspiration probe, complete of 180 cm long connecting tube with female connection. Patient aspiration using only one hand. Rigid canula instantly ready and easy to use. Vacuum regulation by the thumb of the hand. Suitable for every kind of medical secretion suction unit. Singularly sterile packaged.

126140107191 - RIGID PROBE COMPLETE
126140105C01 - YANKAUER
126140104021 - CONNECTING TUBE



SUCTION UNIT CATHETERS

Transparent suction catheter, with silicone. 2 lateral holes, disposable in sterile package with colored connector, available in different sizes.

11214101003 - CATHETER BLACK CH10
11214101004 - CATHETER WHITE CH12
11214101005 - CATHETER GREEN CH14
11214101006 - CATHETER ORANGE CH16
11214101007 - CATHETER RED CH18
11214101008 - CATHETER YELLOW CH20



CATHETER JOINT FINGER TYP - 126140108001

Suction-catheter joint "FINGER TYP" supplied in sterile package. Hole for vacuum regulation.



SILICONE PATIENT TUBE - BSU902

Autoclavable medical silicone patient tube
Length: about 130 cm

ACCESSORIES FOR JARS FOR BOSCAROL SUCTION UNIT



SERRES DISPOSABLE BAG - 57157

SERRES disposable cartridge 1000 ml with lid and integrated antibacterial filter. Exclusively compatible with suction units equipped with OB-J bottle (BSU506). Class I EC medical device.



PROTECTION FILTER FOR OB-J FA JAR

Disposable protection filter with PTFE filter holder and hermetic polypropylene container.

Max. applicable pressure up to 1 bar (100kPa). Mandatory for suction units as per ISO 10079-1:1999.

It is a CE medical device.

Device sold individually and available in two codes BSU999 and M03.1.003. The two codes have similar technical features and functions.

BSU999 - PROTECTION FILTER
M03.1.003 - PROTECTION FILTER

OTHERS SUCTION UNITS



STERILE NEWBORN SUCTION UNIT - 400100

Sterile newborn suction unit.
Lightweight, economic and efficient.
Suitable inside newborn reanimation cases
and sets. Disposable.

RECHARGING AND SUPPLYING SYSTEM FOR SUCTION UNITS



OB WB WALL BRACKET - BSU810

Wall bracket anchorage for OB2012, OB3000 and OB1000 medical suction units. Anchorage system which integrates the recharging system. Anchorage system tested up to 10g capacity as obligatory as per EN 1789 directive. Complete with cable for the connection to the ambulance electric system. 12V power supply for OB2012 and OB1000; 11÷30V power supply for OB3000. Weight: about 0,8 kg
Dimensions: about 240 x 255 x 5 mm



MAIN POWER SUPPLY ADAPTER 2 P 14,0 V

2-pole male connection cable complete with retaining ring for BSU. Wide range input voltage: 100÷240VAc can be used at 50 or 60 Hz. Output voltage 14 Vcc. Delivered rated power: 5 Adc. It can be connected the suction unit while operating on the patient. It is an accessory of medical device and is CE marked. Available in EU, UK and JP version.

BSU895EU - POWER SUPPLY 14,0 V 2P - EU
BSU895UK - POWER SUPPLY 14,0 V 2P - UK
BSU895JP - POWER SUPPLY 14,0 V 2P - JP



RECHARGING-SUPPLYING CABLE 2P 12 V - BSU855

2 poles recharging and power supply cable for Medical Suction Units models. Supplied with universal car lighter plug.



PRACTICAL AND RESISTANT


Our rescue rucksacks are designed to guarantee practicality and resistance. We personally check the quality of our productions to provide even better products.

EMERGENCY RUCKSACKS



EMERGENCY RUCKSACK VIPER 2.0 - BOR87000N

Red tearproof and water-proof material. Double handle (upper or lateral) for easier transport. Large front opening pocket with 2 internal elastic bands. 2 small side pockets and a large one on the outside. Front, side and rear reflective strips (on padded shoulder straps). The part in contact with the back of the user is equipped with cushions. 8 removable coloured pouches. Internal central partition with zipper, complete with elastic bands. Reinforced bottom and rubber feet. Dimensions: about 46x26x50cm. Weight: about 3.8 kg.

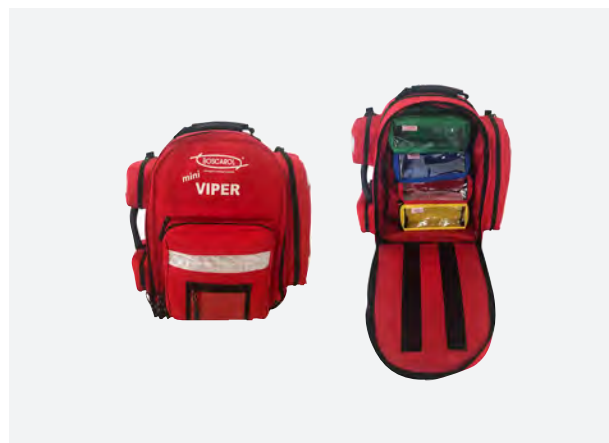

EMERGENCY RUCKSACK MINI VIPER 2.0 - BOR87002N

Red tearproof and water-proof material. Double handle (upper or lateral) for easier transport. Large front opening pocket with 2 internal elastic bands. 2 small side pockets and a large one on the outside. Front, side and rear reflective strips (on padded shoulder straps). The part in contact with the back of the user is equipped with cushions. 4 removable coloured pouches. Reinforced bottom and rubber feet. Dimensions: about 44x26x44 cm. Weight: about 2.3 kg.



EMERGENCY RUCKSACK MASTERMED 2.0 - BOR87020N

Red tearproof and water-proof material. Double handle (upper or lateral) for easier transport. Two front pockets, 1 equipped with elastic bands, made with transparent internal compartments. Front, side and rear reflective strips (on padded shoulder straps). The part in contact with the back of the user is equipped with cushions. 8 removable coloured pouches. Dimensions: about 35x30x53 cm. Volume: about 43 litres. Weight: about 4.3 kg.



EMERGENCY RUCKSACK MAXIMED - BOR87010N

Red tearproof and water-proof material. Front, side and rear reflective strips (on padded shoulder straps). The part in contact with the back of the user is equipped with cushions. 2 front pockets, 2 side pockets, one with 4 coloured pouches. 4 removable coloured pouches. Oxygen cylinder holder straps placed centrally. Dimensions: 48x30x52 cm. Volume: about 54 litres. Weight: about 4,5 kg.

EMERGENCY BAGS



33 BAG - BOR87300

Made of red tearproof and impermeable fabric material. 1 capacious outer pocket. 5 inner adjustable spaces. Inside provided with elastics and a documents pocket. 1 transparent frontal pocket for association logo.

Dimensions: about 36x18x35(h) cm.

Weight: about 1,6 kg.



PROFI BAG MINI - BOR89000

Made of impermeable and tearproof red fabric material. With handle and shoulder straps. Modular interior. 2 wide lateral external pockets. Internal pocket with transparent.

Dimensions: about 36x26x20(h) cm.

Weight: 0,7 kg.



PROFI BAG - BOR89010

Made of impermeable and tearproof red fabric material. With handle and shoulder straps. Modular interior. 1 frontal pocket. Arranged to hold 1-2 lts oxygen cylinder. 2 wide lateral external pockets. Internal pocket with transparent.


Dimensions: about 60x25x27(h) cm.

Weight: 1,4 kg.

TRAUMA BAG


TRAUMA BAG 10 - BOR89040

Easy to carry either by hand, on shoulder or as a backpack
4 external lateral pockets, 2 with elastics, 4 inner modular spaces,
1 inner removable pouch with transparent frontal
Dimensions: about 45x33x30(h) cm.
Weight: about 2,2 kg.

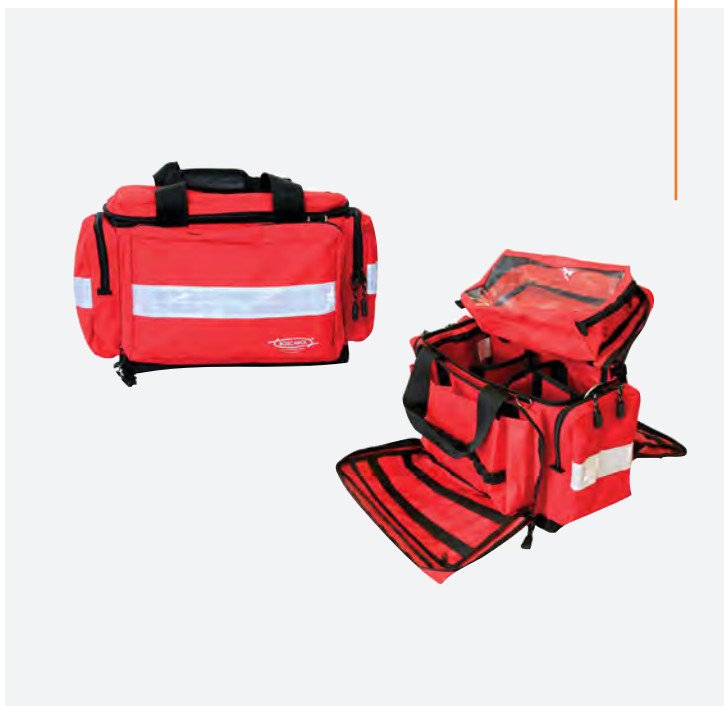

TRAUMA BAG 20 - BOR89050

Easy to carry either by hand, on shoulder or as a backpack. 4 external lateral pockets, 2 with elastics, 5 inner modular spaces, 1 inner removable pouch with transparent frontal.
Dimensions: about 53x33x34(h) cm
Weight: about 2,2 kg

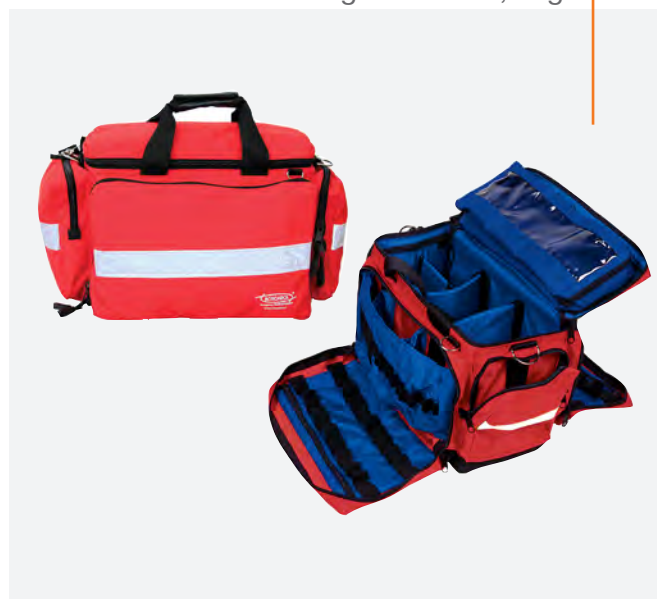

TRAUMA BAG 30 - BOR89060

Easy to carry either by hand, on shoulder or as a backpack. 5 external lateral pockets, 3 with elastics, 5 inner modular spaces, 2 inner removable pouches with transparent frontal.
Dimensions: about 63x37x31(h) cm
Weight: about 2,3 kg.



BOR89040 - TRAUMA BAG 10 RED
BOR89050 - TRAUMA BAG 20 RED
BOR89060 - TRAUMA BAG 30 RED

Pediatric edition of TRAUMA BAG is available (p. 29)

EMERGENCY RUCKSACK ATLS 2.0

ATLS - BOR81145N


Inner, rigid basin with central oxygen cylinder holder and 4 internal removable coloured pouches. The part in contact with the back of the user is equipped with cushions.

Dimensions of all ATLS models:

about 32x25x53cm



Volume: about 35 litres.

Weight: about 3.3 kg.

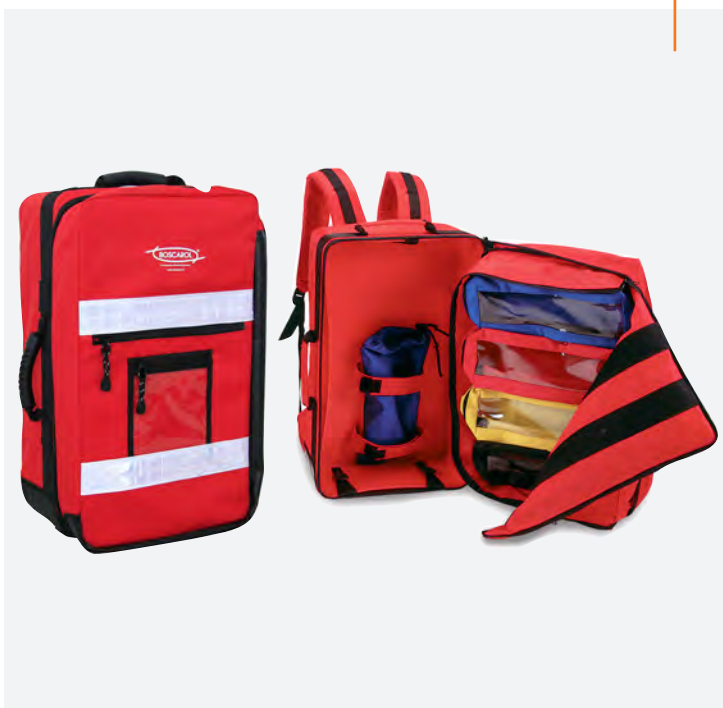

ATLS MODULE - BOR81146N

Inner, rigid basin with 4 compartments and 4 coloured removable pouches. Internal central partition with zipper, and 4 internal removable coloured pouches. The part in contact with the back of the user is equipped with cushions.

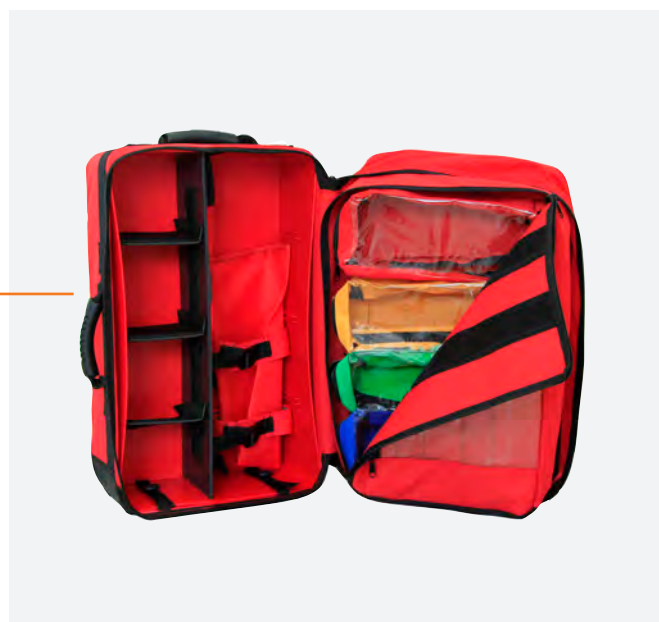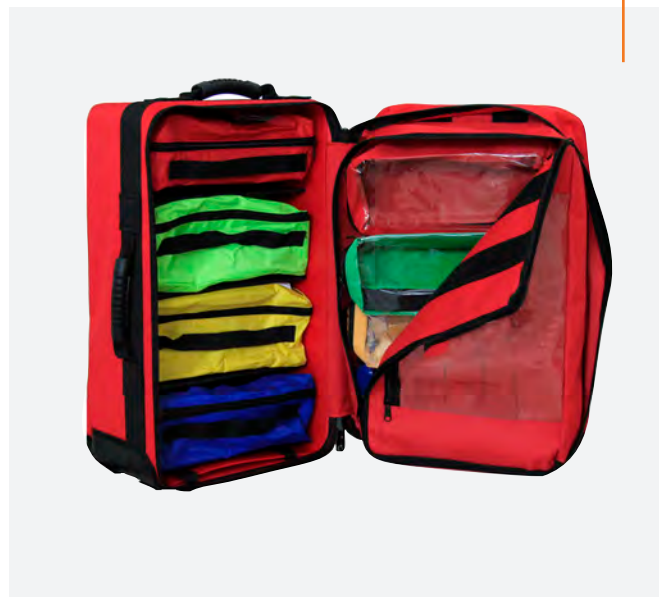

ATLS WK - BOR81147N

Inner, rigid basin with 4 compartments (removable partitioning modules) and 1 lateral oxygen cylinder holder. Internal central partition with zipper, and 4 internal removable coloured pouches. The part in contact with the back of the user is equipped with cushions.

BOR81145N - EMERGENCY RUCKSACK ATLS 2.0
BOR81146N - EMERGENCY RUCKSACK ATLS MODULE 2.0
BOR81147N - EMERGENCY RUCKSACK ATLS WK 2.0

EMERGENCY BAGS



MINI RESCUE BAG - BOR80290

Compact and practical emergency bag. Hand and shoulder transportable. Made in tearproof and waterproof material. 2 windows for the insertion of identity cards/logos. Lower part with inner shockproof basin. Upper part with elastics. Size: about 39 x 18 x 25h cm. Weight: about 1 kg.



RESCUE BAG MODULE - BOR80340

Compact and practical empty rescue bag. Impermeable tearproof material. Easy to carry either by hand, on shoulder or as a backpack. 2 frontal transparent pockets for association logo. Bottom side with shockproof inner basin. 8 removable pouches with transparent. Upper part with inner elastics suitable for equipment fixing. Dimensions: about 47x25x27(h) cm. Weight: about 2,4 kg.



TRAUMA BAG 10 PEDIATRIC - BOR89040P

Made of tear-proof and waterproof nylon. The fabric has bright colours and textures suitable for use in the paediatric sector. Completely removable and washable according to the specified conditions. 4 reflective strips placed on the sides of the bag. Handles, belt and shoulder straps for comfortable and safe transport. 2 external side pockets with internal spaces and elastic bands. The inside of the bag is divided into 4 separate spaces that can be removed. In the upper part there is a removable bag with a transparent front. Dimensions: about 44x30x29 cm. Weight: about 2.1 kg.

DRUGS BAGS



DRUGS BAG 28 PHIALS - BOR89330

Book-opening drugs case empty. Made in tearproof impermeable red fabric material. 2 rear Velcro-strips for the fixing in bags and emergency rucksacks. 2 removable plates with elastics suitable for 28 different sized phials. Central padding for ampoules record. Dimensions: about 25x12x7(h) cm. CE medical device.



DRUGS BAG 84 PHIALS - BOR89332


Big book-opening drugs bag. 2 parts with adjustable elastics to feature 84 different sized phials. Central padding with transparent pocket for ampoules list and expiry date. Made in red, tear-resistant, water-resistant material. 2 external Velcro-strips for the fixing in bags or rucksacks. Dimensions: about 25x25x7h cm. CE medical device.

POUCHES



BLUE DOUBLE POUCH - BOR81042

Dimensions: about 25x25x8 cm

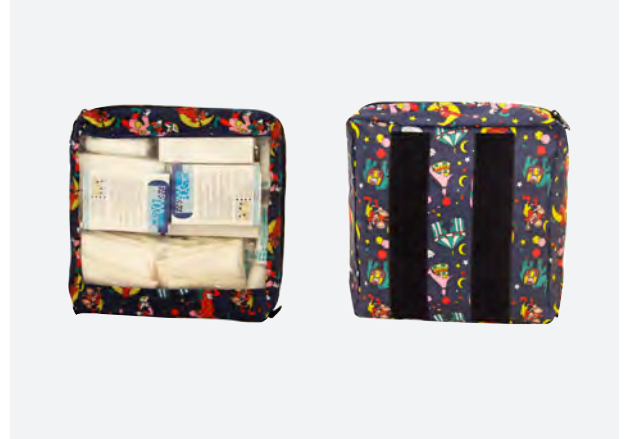

FANTASY POUCH - BOR81056

Dimensions: about 25x25x8 cm



RED PVC MINI POUCH - BOR89469

Dimensions: about 10x7x13 cm
Designed to be attached to the OB2012,
OB1000 AND OB MINI Boscarol Suction Unit



STANDARD POUCHES

Dimensions: about 27x12x9 cm

BOR81024 - STD POUCH PINK
BOR81028 - STD POUCH RED
BOR81030 - STD POUCH GREEN
BOR81032 - STD POUCH YELLOW
BOR81034 - STD POUCH BLUE
BOR81036 - STD POUCH PURPLE

POUCHES



BIG POUCHES WITH HANDLE

Dimensions: about 30X12X13

BOR89440 - GREEN BIG POUCH WITH HANDLE
BOR89442 - RED BIG POUCH WITH HANDLE
BOR89444 - BLUE BIG POUCH WITH HANDLE
BOR89446 - YELLOW BIG POUCH WITH HANDLE



POUCHES WITH HANDLE

Dimensions: about 23x14x11

BOR89420 - YELLOW POUCH WITH HANDLE
BOR89422 - BLUE POUCH WITH HANDLE
BOR89424 - GREEN POUCH WITH HANDLE
BOR89426 - RED POUCH WITH HANDLE



VIPER POUCHES

Dimensions: about 19x12x6 cm
Velcro-strips: female

BOR81058- YELLOW VIPER POUCH
BOR81060- GREEN VIPER POUCH
BOR81062- RED VIPER POUCH
BOR81066- VIOLET VIPER POUCH
BOR81068- BLACK VIPER POUCH
BOR81070- ORANGE VIPER POUCH



MINI VIPER POUCHES

Dimensions: about 19x12x6 cm
Velcro-strips: male

BOR82000 - VIOLET MINI VIPER POUCH
BOR82001 - ORANGE MINI VIPER POUCH
BOR82002 - BLACK MINI VIPER POUCH
BOR82006 - DARK BLUE MINI VIPER POUCH

POUCHES



MINI POUCHES

Dimensions: about 15x12x8,5 cm

BOR89458 - MINI POUCH BLACK
BOR89460 - MINI POUCH RED
BOR89462 - MINI POUCH BLUE
BOR89464 - MINI POUCH GREEN



FABRIC BAG - BOR84072

Fabric empty bag. Tearproof impermeable blue material. With shoulder strap.
Dimensions: about 29x16x16 cm.
Weight: about 210 g. Supplied empty.



ROLL-CASE FOR LARYNGOSCHOPE - BOR84041

Roll-case for laryngoscope - red
5 inner compartments.
Made in washable cloth with inner stuffing.
Dimensions: about 19x34 cm.

EMPTY/COMPLETE E-HOLSTER AND HOLSTER



EMPTY E-HOLSTER

E-Holster belt made in tearproof and water-resistant red material. 5 compartments with pockets for gloves and for Pen Light. Belt loop of about 5,5 cm. Dimensions: about 14x3x19(h) cm. Weight: about 100 g. Supplied empty.

BOR84062 - BLUE EMPTY E-HOLSTER
BOR84063 - RED EMPTY E-HOLSTER



COMPLETE E-HOLSTER

E-Holster belt made in tearproof and water-resistant red material. 5 compartments with pockets for gloves and for Pen Light. Belt loop of about 5,5 cm. Complete with: 1 first aid scissors, 1 anatomic forceps, Cocher forceps, 1 Pen Light, 1 Tourniquet. Dimensions: about 14x3x19(h) cm. Weight: about 200 g.

BOR84064 - BLUE COMPLETE E-HOLSTER
BOR84065 - RED COMPLETE E-HOLSTER



HOLSTER PROFESSIONAL - BOR206

Particularly functional belt case Made in tearproof impermeable fabric material. Suitable to contain one sphygmomanometer. 3 inner spaces with closing. Dimensions: about 34x11x8(h) cm. Weight: about 100 g. Supplied empty.

CABINETS



VISIO DEFIBRILLATOR METAL CABINET - DEF040

Metal Cabinet for Defibrillator with see through door. Made of Galvanized Carbon Steel PO2. Enameled by Epoxy Polyester white powder RAL9010 Heated at 180°. Transparent door with opening at 180° with a galvanized strong pivot. Strong Nickel plated zinc alloy closing device. Can be sealed or locked. Back holes for wall hanging. Dimensions: 425x425x160 mm. Weight: 6 kg. Supplied empty.



CABINET - CAV605

Carbon steel cabinet, coated in epoxy polyester. Rounded edges. Three steel shelves. Door with 3 accessories holder. Locker with key. Dimensions: 402x202x483(h) mm. Weight: 6,60 kg. Supplied empty.

ALUMINIUM CASES



ALUMINIUM CASE TYPE A - VAL42000

Type A emergency medical case made of anodized aluminium. High resistance and long durability. Inside made in washable and disinfectable thermo-format with modular system. Arranged to accommodate oxygen cylinders (1 or 2 litres).

Dimensions: about 52x36x17(h) cm.

Weight: about 5,7 kg. Supplied empty.

VAL42000 - ALUMINIUM CASE TYPE A
VAL42010 - 2 AMPOULE HOLDER RAIL
DEPOT FOR CASE TYPE A



ALUMINIUM CASE TYPE B - VAL43000

Type B emergency medical case, made of anodized aluminium. High resistance and long durability. Inside made in washable and disinfectable thermo-format with modular system. Arranged to accommodate oxygen cylinders (1 or 2 litres). Removable internal basin made in thermo-format with rail ampoules holder (max 60 pcs).

Dimensions: about 52x36x24(h) cm.

Weight: about 7,6 kg. Supplied empty.



RAIL AND EASY RAIL VIALS HOLDER

We have available RAIL and EASY-RAIL vials holders made of painted metal with internal rubber, suitable for vials of any size. Available sizes indicated on request.

CASES IN ABS



SAN CASE IN ABS - VAL40130

San type case empty. Made in orange shockproof ABS material. Wall bracket anchorage system. 2 equal sides made of 2 adjustable inner compartments. 2 transparent dividers. Dimensions: about 32,5x12,5x21,5(h) cm. Weight: about 1,1 kg.



MULTI CASE IN ABS - VAL40110

Multi type case empty. Made in orange shockproof ABS material. Wall bracket anchorage system. 2 equal sides made of 2 adjustable inner compartments. 2 transparent dividers. Dimensions: about 42x15x30(h) cm. Weight: about 3,3 kg.



MEDIC CASE IN ABS

Medic 2 case empty with orange or transparent cover. Made in shockproof orange ABS material. Equipped with wall bracketing. Dimensions: about 39,5x13,5x26(h) cm. Weight: about 1,6 kg.

CAV399 - MEDIC 2 CASE TRANSPAR. COVER
CAV401 - MEDIC 2 CASE ORANGE COVER



FLEXI CASE IN ABS - CAV406

Flexi case empty. Made in shockproof orange polypropylene material. Equipped with wall bracketing. Dimensions: about 25,5x9x19,5 (h) cm. Weight: about 0,45 kg.

WATER PROOF CASES



VAL43140 - SMALL BOSCAROL HIGH IMPACT CASE
VAL43141 - MEDIUM BOSCAROL HIGH IMPACT CASE
VAL43142 - LARGE BOSCAROL HIGH IMPACT CASE

BOSCAROL HIGH IMPACT CASE

Boscarol High Impact small orange case with high resistant thick walls. Automatic pressure release valve. Easy to open double throw latches. Load tested soft grip handle. Holes for double padlocking. Holes for shoulder strap. Stackable interlocking system. Thick body material. Extra fixing points. Fully length hinges. Nylon pins, corrosion free. Pick N Pluck foams in the bottom and convoluted foam available on demand. Screen printing space on the lid. Rectangular resin label space on the lid.

Available in 3 sizes:

Small: about 33,6x30x14,8 cm.

Weight: 1,75 kg.

Medium: 46,4x36,6x17,6 cm.

Weight: 2,78 kg.



EXPLORER CASE

Emergency waterproof and watertight case up to 10 meters deep. With pressurization valve. Sturdy closing system with predisposition for a lead lock.

Available in different sizes.

EXPLORER MINI

Dim.: 30,5x19,5x27h cm - Weight: 1,9 kg

EXPLORER P

Dim.: 41x20,5x34 cm - Weight: 3,1 kg

EXPLORER M

Dim.: 52x23x43,5 cm - Weight: 4,6 kg

EXPLORER

Dim.: 44,5x34,5x19 cm - Weight: 4 kg

EXPLORER G

Dim.: 65x24x51 cm - Weight: 7,1 kg

VAL43151 - EXPLORER MINI
VAL43152 - EXPLORER P
VAL43162 - EXPLORER M
VAL43164 - EXPLORER
VAL43172 - EXPLORER G
VAL43174 - EXPLORER P WITH 4 POUCHES
VAL43176 - EXPLORER M WITH 5 POUCHES
VAL43178 - EXPLORER G WITH 6 POUCHES



PRODUCTS FOR HEALTH EMERGENCY

Our immobilization and transport medical devices meet the standards in force when expected and were born from the desire to provide products dedicated to the specific sector of health emergency.

MADE IN ITALY

The MADE IN ITALY manufacturing allows us to offer high quality medical devices, by verifying and tracing every phase of the production process.

SUCTION VALVE



All our vacuum medical devices are equipped with a suction valve specifically designed and patented by Oscar Boscarol srl.

VACUUM MATTRESSES

SDM - IMM120012



SDM vacuum mattress, made of high-strength fabric, intended for the immobilisation and transport of polytrauma patients. Equipped with 8 handles for lifting and transportation. Suction valve located close to the head on the outer side, easily accessible even with immobilized patient. Equipped with 3 immobilisation straps. Radio-transparent. Equipped with 3 internal suction chambers for better adaptation to the anatomical shape of the patient. Dimensions: 200x80x6h cm

Weight: 7 kg

In compliance with the standard EN 1865-1
It is a CE medical device in compliance with the European Regulation MDR 2017/745.

IMM120012 - SDM VACUUM MATTRESS
XIMM120014 - SDM VACUUM MATTRESS
COMPLETE WITH BAG AND PUMP

PROFI - IMM120150



PROFI vacuum mattress, made of high-strength fabric, intended for the immobilisation and transport of polytrauma patients. The mattress includes a semi-rigid support in which 8 handles for lifting and transport are located. Suction valve located close to the head on the outer side, easily accessible even with immobilized patient. Equipped with 3 immobilisation straps. Radio-transparent. Equipped with 3 internal suction chambers for better adaptation to the anatomical shape of the patient. Dimensions: 205x95x10h cm

Weight: 10 kg

In compliance with the standard EN 1865-1
It is a CE medical device in compliance with the European Regulation MDR 2017/745.

IMM120150 - PROFI VACUUM MATTRESS
XIMM120130 - PROFI VACUUM MATTRESS
COMPLETE WITH BAG AND PUMP

VACUUM MATTRESSES

PROFI EVOLUTION - IMM120129



PROFI EVOLUTION vacuum mattress, made of high-strength fabric, intended for the immobilisation and transport of polytrauma patients. The mattress includes a semi-rigid support in which 8 handles for lifting and transport are located. Suction valve located close to the head on the outer side, protected by the semi-rigid support, easily accessible even with immobilized patient. Equipped with two specific holes for fixing the mattress to metal stretchers used in air rescue operations. Equipped with 3 immobilisation straps. Radio-transparent. Equipped with 3 internal suction chambers for better adaptation to the anatomical shape of the patient.

Dimensions: 205x95x10h cm

Weight: 8.2 kg

In compliance with the standard EN 1865-1
It is a CE medical device in compliance with the European Regulation MDR 2017/745.

IMM120129 - PROFI EVOL. VACUUM MATTRESS
XIMM120131 - PROFI EVOL. VACUUM
MATTRESS COMPLETE WITH BAG AND PUMP



BABY - IMM120060

BABY vacuum mattress, made of high-strength fabric, intended for the immobilisation and transport of children. Equipped with 4 handles for lifting and carrying. Suction valve located close to the head on the outer side, easily accessible even with immobilized patient. Radio-transparent.

Dimensions: 87x65x4h cm

Weight: 1.6 kg

Max. load capacity: 18 kg

VACUUM SPLINT

VACUUM SPLINT



High resistant and washable material. Quick and painless immobilization of all kind of fractures (even exposed ones). Immobilizes without constrictions thus allowing blood circulation. External aspiration valve. 100% X-ray transparent (any metal parts included). Suitable for use on baby patients.

CE medical device.

Dimension hand vacuum splint: 50,5x33 cm.

Dimension arm vacuum splint: 75x47 cm.

Dimension leg vacuum splint: 92x64 cm.

Dimension leg 108 vacuum splint: 108x64 cm. Singularly or in set packaged, complete with transport bag or transport bag and suction pump.

IMM120172 - VACUUM SPLINT LEG 108CM

IMM120174 - VACUUM SPLINT HAND

IMM120176 - VACUUM SPLINT ARM

IMM120178 - VACUUM SPLINT LEG

IMM120200 - SET 2 VACUUM SPLINT WITH TRANSPORT BAG

IMM120180 - SET 2 VACUUM SPLINT COMPLETE (TRANSPORT BAG AND SUCTION PUMP)

IMM120190 - SET 3 VACUUM SPLINT WITH TRANSPORT BAG

IMM120184 - SET 3 VACUUM SPLINT COMPLETE (TRANSPORT BAG AND SUCTION PUMP)

IMM120218 - TRANSPORT BAG FOR VACUUM SPLINT

ACCESSORIES FOR VACUUM DEVICES



ALUMINIUM MAXI SUCTION PUMP - IMM120240

Suction pump for vacuum mattress and splints
Made in aluminium with foot support for easy aspiration.

Dimensions: about 13x5,5x54(h) cm

Weight: about 1,4 kg



ALUMINIUM MINI SUCTION PUMP - IMM120260

Suction pump for vacuum mattress and splints. Made in aluminium.

Dimensions: about 4x32 cm

Weight: about 500 g



TRANSPORT BAG FOR VACUUM MATTRESS - IMM120030

Transport bag for vacuum mattress.
Tearproof impermeable material.
Transport handles.

Dimensions: about 89x24x72(h) cm



REPAIR KIT- AIM5000

It contains five pieces of pvc fabric and a bottle of special glue: make sure the damaged spot is not covered with oil or grease, clean it thoroughly and make it rough. Apply the special glue and let it dry (check with the back of the fingers). Apply the fabric with firm pressure to that area. Check the seal.

RIGID SPLINT AND SSPLINT

RIGID SPLINT

Rigid splint for upper/inferior limbs immobilization. Made in PVC and Nylon. Washable and X-ray transparent. Velcro closing. Dimensions upper limbs: about 51x27,5x3,5 cm. Dimensions inferior limbs: 74x40x3,5 cm. Available in set of 2 vacuum splints complete with transport bag.

IMM123100 - UPPER LIMBS RIGID SPLINT
IMM123110 - INFERIOR LIMBS RIGID SPLINT
IMMM123120 - SET 2 RIGID SPLINT WITH BAG
IMM123125 - TRANSPORT BAG FOR RIGID SPLINT


SSPLINT

Splint system set with aluminium core for easy limbs immobilization in case of fracture or sprain. Made in neoprene material. Adaptable to fit each body frame. It does not cause tissue compressions. It maintains body temperature. It reduces congestion around injured areas. Velcro modelling closing. Lightweight and easily washable even in washing-machine. Least encumbrance. Dim. leg ssplint: 58x21x0,6 cm (not suitable for femur fractures). Dim. arm ssplint: 60,5x17x0,6 cm. Dim. forearm ssplint: 39x16x0,6 cm. Dim. wrist ssplint: 32x15x0,6 cm. Dim. ankle ssplint: 57x18x0,6 cm. Set of 5 vacuum splints complete with bag.

IMM124030 - SSPLINT OB SET 5 PCS WITH BAG
IMM124031 - SSPLINT OB LEG
IMM124032 - SSPLINT OB ARM
IMM124033 - SSPLINT OB FOREARM
IMM124034 - SSPLINT OB WRIST
IMM124035 - SSPLINT OB ANKLE
IMM123154R - RED BAG FOR SSPLINT

ADAPTABLE AND BREATHABLE FINGER SPLINT - 50052

Suitable for small bone detachments, finger dislocations, phalanx fractures, contusive traumas. Class I non-sterile medical device in compliance with Directive 93/42/EC and following amendments and integrations. Size: M/L.

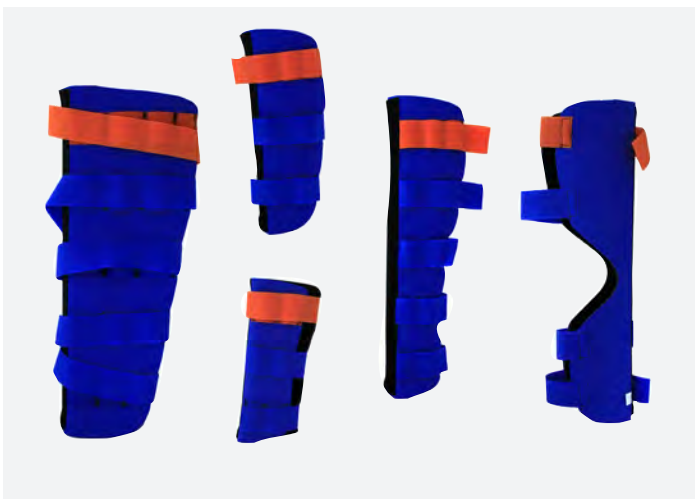

OB NEW ONE-PIECE CERVICAL COLLARS

One-piece cervical collar, made with injection mould system. Perfect immobilization of cervical section. Frontal opening for trachea check. Made of polypropylene X-ray transparent with non absorbent padding. CE medical device. Available in different sizes and in set of 2,3, 4 or 6 collars complete with transport bag.




0 INFANT PINK - IMM124100

Size: 0 - Infant pink



1 PEDIAT. LIGHT BLUE - IMM124102

Size: 1 - Pediatric light blue

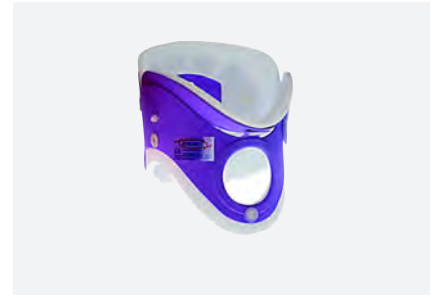

2 NECKLESS VIOLET - IMM124104

Size: 2 - Neckless violet



3 SHORT BLUE - IMM124106

Size: 3 - Short blue



4 REGULAR ORANGE - IMM124108

Size: 4 - Regular orange



5 TALL GREEN - IMM124110

Size: 5 - Tall green



SET OF 3 OR 4 COLLARS

IMM124112: set of 3 OB NEW one-piece collars (size from 1 to 3) with transport bag.

IMM124114: set of 4 OB NEW one-piece (size from 1 to 4) with transport bag.



SET OF 6 COLLARS - IMM124116

Set of 6 OB NEW one-piece collars (sizes from 0 to 5), with transport bag.

IMM120656: transport bag.



SET OF 2 COLLARS - IMM124118

Set di 2 OB NEW one-piece collars (sizes 0 and 1), with transport bag.

OBNECK TWO-PIECES CERVICAL COLLARS

Two-pieces cervical collar, made with injection mould system. Perfect immobilization of cervical section. Frontal opening for trachea check. Made of polypropylene X-ray transparent with non absorbent padding. CE medical device. Available in different sizes and in set of 2,3 or 6 collars complete with transport bag.



0 INFANT PINK - IMM120782

Size: 0 - Infant pink



1 PEDIAT. LIGHT BLUE - IMM120784

Size: 1 - Pediatric light blue



2 NECKLESS VIOLET - IMM120786

Size: 2 - Neckless violet



3 SHORT BLUE - IMM120788

Size: 3 – Short blue



4 REGULAR ORANGE - IMM120790

Size: 4 - Regular orange



5 TALL GREEN - IMM120792

Size: 5 - Tall green



SET 3 COLLARS - IMM120794

Set of 3 OBNECK two-pieces collars (sizes from 1 to 3), complete with transport bag.



SET 6 COLLARS - IMM120796

Set of 6 OBNECK two-pieces collars (sizes from 0 to 5), complete with transport bag.



SET 2 COLLARS - IMM120798

Set of 2 OBNECK two-pieces collars (sizes 0 and 1), complete with transport bag.

IMM120653 - TWO-PIECES COLLARS TRANSPORT BAG
IMM120658 - FANTASY TWO-PIECES COLLARS TRANSPORT BAG

ADJUSTABLE CERVICAL COLLARS

ADJUSTABLE COLLARS CHAMPION



Double valve collar adjustable in different heights. Allows immobilization of the cervical spine, minimizes flexion-extension and ensures better control of axial rotation and lateral flexion movements. Central hole for the control of the trachea and the execution of a tracheotomy procedure. The end of the Velcro is in orange color for a better identification of the grip. Polyethylene structure, black LD padding, nylon buttons and Velcro strap latch. Radiotransparent - Latex-Free - Unsterile - Single-use - Disposable.

9030 - CHAMPION COLLAR ADULT

9031 - CHAMPION COLLAR PEDIATRIC

ADJUSTABLE COLLARS FAST



Cervical collar adjustable to different heights. Allows immobilization of the cervical spine and minimizes flexion-extension. Central hole for the control of the trachea and the execution of a tracheotomy procedure. Polyethylene structure, black LD padding, nylon buttons and Velcro strap latch. Radiotransparent - Latex-Free - Unsterile - Single-use - Disposable.

9015 - FAST COLLAR ADULT

9015PED - FAST COLLAR PEDIATRIC

ADJUSTABLE COLLARS PATRIOT



Cervical collar 4 positions adjustable (Low, Short, Regular and Tall). Simple and versatile usage. Excellent immobilization. Wide frontal opening allowing an adequate access to windpipe. Easily washable. X-Ray transparent.

PHP-PA - PATRIOT ADULT

PHP-PP - PATRIOT PEDIATRIC

PHP-PA-G - PATRIOT ADULT OLIVE GREEN

CERVICAL COLLARS AND IMMOBILIZATION DEVICES



STIFNECK COLLARS

Stifneck cervical collar with 4 sizes (Low, Short, Regular e Tall). Frontal opening for windpipe check. Easy to wash. X-ray transparent.

IMM120648 - STIFNECK ADULT COLLAR
IMM120650 - STIFNECK PEDIATRIC COLLAR



ECO COLLARS

Cervical collar made in high density Elastolen/Veolene, non-deformable (closed cells polyethylene) resistant to external agents, totally non-toxic. Velcro closing. X-ray transparent.

IMM120776 - ECO COLLAR SMALL N. 0
IMM120774 - ECO COLLAR MEDIUM N. 2
IMM120778 - ECO COLLAR BIG N. 3
IMM120770 - SET 3 ECO COLLARS (0-2-3)



SOFT COLLARS

Soft collars for cervical immobilization. Structure in high density foam rubber with tubular cotton jersey covering. Velcro closing.

IMM120830 - SOFT COLLAR SIZE S
IMM120832 - SOFT COLLAR SIZE M
IMM120834 - SOFT COLLAR SIZE L
IMM120836 - SOFT COLLAR SIZE XL



RESCUE CAN - 6400300

Rescue can, rescue buoy, light and imperishable. Printed in high density polyethylene without joints. Endowed with floating 186 cm rope and shoulder strap. Non-homologated by nautical authorities. It is not a medical device.

Dimensions: about 23x16x68 cm.

Weight: about 1,2 kg.

OB SPINAL BOARD AND TRANSFER SHEET



OB SPINAL BOARD

High resistant plastic spinal board, suitable on difficult emergency activities. 4 buttonholes for quick and total fixing of head immobilizer and 2 hollows placed on leaning side for a traditional base blocking thus allowing strap belts not to get teared during use. Provided with 20 transport handles. X-ray translucent.

Maximum load: about 150 kg

Dimensions: about 184x41x5,5 cm

Weight: about 5,5 kg

IMM121074 - OB SPINAL BOARD

IMM121076 - OB SPINAL BOARD WITH SPIDER BELTS AND OB HEAD IMMOBILIZER

IMM121075 - OB SPINAL BOARD WITH SPIDER BELTS AND WIND HEAD IMMOBILIZER



TRANSFER SHEET - IMM120321

New transfer sheet, even more resistant. Made of red tearproof impermeable PVC material. 8 side handles. Easily washable and disinfectable. Perimeter reinforced with a 4cm, double sewn band on both sides.

Dimensions: about 187x70 cm

Weight: about 1200 g

Loading: about 150 kg

In accordance to EN1865-1 norm

UNIVERSAL HEAD IMMOBILIZER



OB HEAD IMMOBILIZER - IMM121628

Made in high density treated yellow foam material, water repellent. Featuring 1 basic pillow, two head support pillows with opening for ears access, 2 chin supports with central padding and 2 straps with Velcro already integrated to the basis for the positioning of the head Immobilizer on the stretcher. X-Ray translucent. Dimensions: about 42,5x26,5x16(h) cm Weight: about 900 g

IMM121628 - OB HEAD IMMOBILIZER
AIM4190 - BASE FOR HOB EAD IMMOBILIZER
AIM4191 - SET OF 2 CHIN STRAPS FOR OB HEAD IMMOBILIZER
AIM4192 - SET OF 2 HEAD-BLOCK FOR OB HEAD IMMOBILIZER




WIND HEAD IMMOBILIZER - IMM121640

Realized in stuffed MDF and covered with washable, disinfectable, water repellent and thermo-welded material. 2 head-fixing-pillows with opening for access to the ears are integrated at the base and are part of the base itself if not used. Allows fixed head immobilizer storage with stretcher. Complete of 2 chin straps and 2 Velcro straps, 8 belts with loop and 2 velcro bands for the fixing to the spinal board. X-ray translucent. Dimensions: about 44,5x30x2,5 cm Weight: about 1,2 kg

IMM121640 - WIND HEAD IMMOBILIZER
IMM121634 - SET OF 2 CHIN-STRAPS FOR WIND HEAD IMMOBILIZER
IMM121652 - SET OF 4 HEAD IMMOBILIZER BELTS
IMM121654 - SET OF 2 VELCRO STRAPS FOR WIND HEAD IMMOBILIZERS

IMMOBILIZATION BELTS

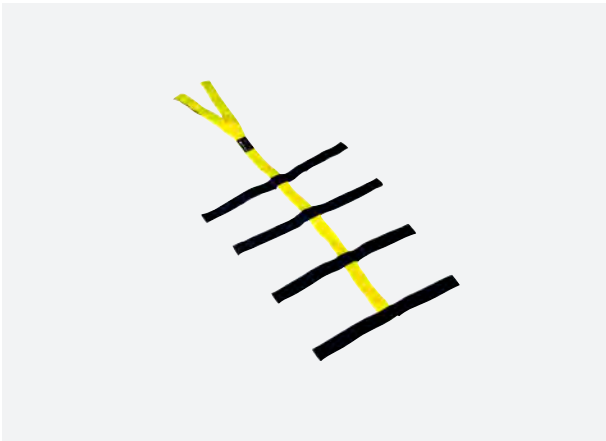

SPIDER IMMOBILIZATION BELTS - IMM121033

Allows the immobilization of an adult or paediatric patient on a generic spine board. Made with 1.2 mm thick polyester webbing to ensure high fastening strength and durability. Thanks to the possibility of moving the straps according to the type of patient, (along the vertical axis of the fastening system) it is possible to secure the patient in a stable and safe way. The straps are fastened to the patient's body thanks to 5 cm wide Velcro strips which make both closing and opening operations easy and safe. Complete with a practical case.



SET 3 PLASTIC BUCKLE FIXING BELTS - IMM120280

Set of 3 stretcher fixing belts. Quick hook-unhook plastic buckle. Complete with case with Velcro closing.



SET 2 METAL BUCKLE FIXING BELTS - IMM120282

Set of 2 closed safety belts
Quick hook-unhook metal buckle
Dimensions: about 5x180 cm




2 SET 2 PLASTIC BUCKLE SAFETY BELTS - IMM120284

Set of 2 2-pieces safety belts with loop
Quick hook-unhook plastic buckle for the use on automatic loading stretcher.
Dimensions: about 5x176 cm

CTS - CHILDREN TRANSPORT UNIT - IMM120050

PEDIATRIC EMERGENCY TRANSPORT

Children transport system with pediatric fabric. Complete of safety-belts allowing to block the newborn in 3 positions, avoiding children usual dangerous and unexpected movements.

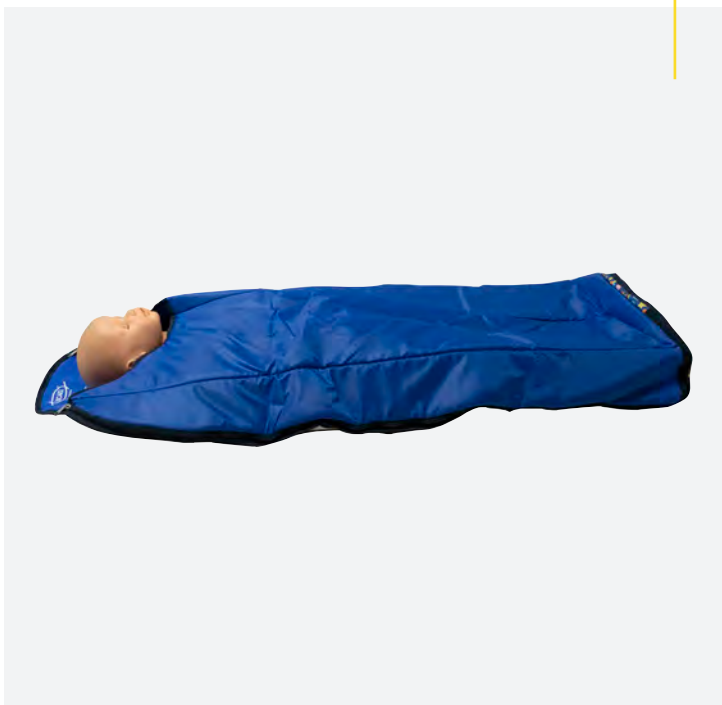

FROM 1 TO 24 MONTHS

The device is intended to be used for new-born and children transport, approximately from 1 to 24 months old (shoulder / crotch length from 35 to 57 cm).



ANCHORAGE TO THE STRETCHER

3 external safety belts allowing child further fastening. 3 belts allowing the fixing of the CTS to the stretcher. Complete of transport bag. CE medical device. Dimensions: about 110x70 cm



INFANT TRANSPORT UNIT - IMM120071

NEWBORN EMERGENCY TRANSPORT

Newborn emergency transportation system, complete with inner fixing belts. Grants a warm, protected and safe environment. 3 transparent sides allowing complete and continuous monitoring of newborn by medical personnel.



INNER FIXING BELTS


Nuova versione PLUS completa di cinghie interne per il fissaggio del paziente. Materasso interno per un trasporto confortevole. Realizzazione a borsa ripiegabile che permette di essere agevolmente richiusa e trasportata.



ANCHORAGE TO THE STRETCHER

2 external handles for transportation. Arranged for the anchorage to the ambulance stretcher. Dimensions open: cm 90x38x34,5(h)
Dimensions closed: cm 90x38x15(h)
Supplied with transport bag.

RESCUE STRETCHER

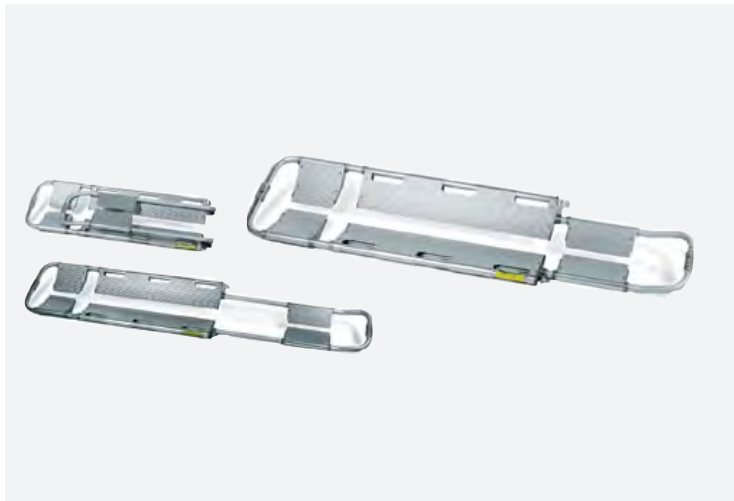

SCOOP STRETCHER - 630

Aluminium scoop stretcher. X-ray translucent. Effective hook/unhook device for stretcher parting. Adjustable length with automatic block. Bag with n. 3 quick unhook safety belts for patient immobilization. Capacity: about 170 kg. Weight: about 10 kg. Height: about 7 cm. Width: about 43 cm. Minimum Length: about 167 cm.

630 - SCOOP STRETCHER
614 - HORIZONTAL STIRRUPS
616 - VERTICAL STIRRUPS



RESCUE STRETCHER TOBOGA - 16100

Rescue stretcher, extremely lightweight, suitable for patient transport in uncomfortable environment conditions. Made in polymeric ABS material with aluminium frame. Complete of expanded polyutherane mattress and adjustable footrest. 4 quick unhook security belts. 12 transport handles. Arranged for the anchorage to helicopters. Dimensions: about 216x61x17h cm. Weight: 15 kg. Capacity: max 300 kg.

16100 - RESCUE STRETCHER TOBOGA
16100-005 - LIFTING HARNESS
16100-002 - FLOTATING SYSTEM

EMERGENCY STRETCHER



EMERGENCY STRETCHER X1

Emergency stretcher made of lightweight aluminium. Foldable in length with safety blocks on use position. Made in resistant fireproof, tearproof and washable material with central reinforce. Complete with 2 automatic unhook safety belts. Removable handles in plastic non-sliding material. Available with feet and wheels.

Dimensions (opened): 230x55,6x13,7 cm.

Dimensions (closed): 192x15x14,5 cm.

Weight: about 9 kg.

IMM121532 - X1 WITH 4 FEETS

IMM121534 - X1 WITH 4 WHEELS

IMM121530 - X1 WITH 2 FEETS AND 2 WHEELS



EMERGENCY STRETCHER X2

Emergency stretcher made of lightweight aluminium. Foldable in length and width with safety blocks on use position. Made in resistant fireproof, tearproof and washable material with central reinforce. Complete with 2 automatic unhook safety belts. Removable handles in plastic non-sliding material. Available with feet and wheels.

Dimensions (opened): 230x55,6x13,7 cm.

Dimensions (closed): 104x15x14,5 cm.

Weight: about 9 kg.

IMM121538 - X2 WITH 4 FEETS

IMM121540 - X2 WITH 4 WHEELS

IMM121536 - X2 WITH 2 FEETS AND 2 WHEELS




ANATOMIC PILLOW FOR STRETCHER - IMM122050

Anatomic pillow for stretcher. Made in antibacterial, fireproof, washable material. Completely thermo-welded.

Dimensions: about 40x27x7 cm.

Weight: about 300 g.

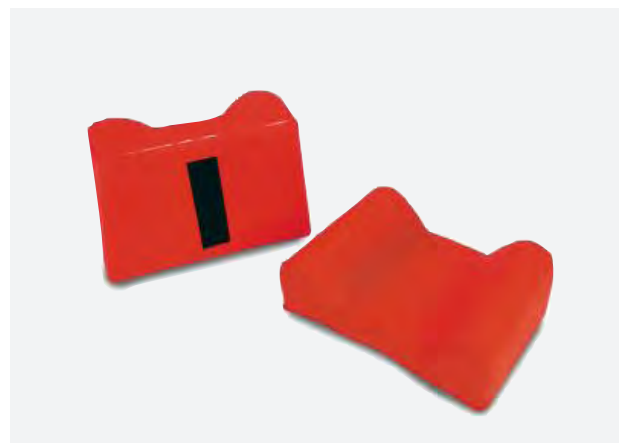

CONFORMING PILLOW FOR STRETCHER - IMM122070

Conforming orange pillow for stretchers. Made in antibacterial, fireproof, washable material. Completely thermo-welded.

Dimensions: about 38x25x10/14h cm.

Weight: about 400 g.



RESUSCITATOR BAGS

SILICONE RESUSCITATOR BAGS



ADULT

Transparent silicone resuscitator bag, autoclavable up to 135°C under the pressure of 2 bar. Volume 1600 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask no. 5, reservoir and oxygen tube Latex free. CE certified in class IIa.

PAL32000 - ADULT SILICONE
PAL32004 - ADULT SILICONE
COMPLETE



PEDIATRIC

Transparent silicone pediatric resuscitator bag, autoclavable up to 135°C under the pressure of 2 bar. Volume 500 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask no. 3, reservoir and oxygen tube. Latex free. CE certified in class IIa.

PAL32010 - PEDIATRIC SILICONE
PAL32014 - PEDIATRIC SILICONE
COMPLETE
PAL32010 - PEDIATRIC SILICONE
ALL-IN-ONE



NEWBORN

Transparent silicone newborn resuscitator bag, autoclavable up to 135°C under the pressure of 2 bar. Volume 240 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask no. 1, reservoir and oxygen tube. Latex free. CE certified in class IIa.

PAL34024 - NEWBORN SILICONE
COMPLETE ALL-IN-ONE

DISPOSABLE RUBBER RESUSCITATOR BAGS



ADULT

Adult non-autoclavable rubber resuscitator bag. Volume 1600 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. To be disinfected in bactericide solution. CE certified in class IIa.

PAL32200 - ADULT RUBBER



PEDIATRIC

Pediatric non-autoclavable rubber resuscitator bag. Volume 500 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask size 3 with reservoir and shockproof oxygen tube. Latex free. To be disinfected in bactericide solution. CE certified in class IIa.

PAL32210 - PEDIATRIC RUBBER
PAL34210 - PEDIATRIC RUBBER
ALL-IN-ONE
PAL34214 - PEDIATRIC RUBBER
COMPLETE ALL-IN-ONE

DISPOSABLE PVC RESUSCITATOR BAGS



ADULT

Disposable economic resuscitator bag. Made in PVC material. Volume 1600 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask size 5 with reservoir and shockproof oxygen tube. Latex free. CE certified in class IIa.

PAL32100 - ADULT PVC
PAL32104 - ADULT PVC COMPLETE



PEDIATRIC

Disposable economic resuscitator bag. Made in PVC material. Volume 500 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask size 3 with reservoir and shockproof oxygen tube. Latex free. CE certified in class IIa.

PAL32110 - PEDIATRIC PVC
PAL32114 - PEDIATRIC PVC
COMPLETE
PAL34110 - PEDIATRIC PVC ALL-IN-ONE



NEWBORN

Disposable economic resuscitator bag. Made in PVC material. Volume 240 ml. Complete with 360° rotating connection patient's valve. Adjustable safety valve to limit the overpressure applied by the patient. Complete version supplied with mask size 1 with reservoir and shockproof oxygen tube. Latex free. CE certified in class IIa.

PAL32120 - NEWBORN PVC
PAL34124 - NEWBORN PVC
COMPLETE ALL-IN-ONE

VENTILATION MASKS



SILICON MASKS

Transparent ventilation mask with silicone autoclavable lip-board ensuring maximum adherence to the face. The transparent shell allows patient's monitoring during ventilation.

PAL32050 - SILICONE MASK SIZE 0
PAL32054 - SILICONE MASK SIZE 1
PAL32058 - SILICONE MASK SIZE 2
PAL32062 - SILICONE MASK SIZE 3
PAL32066 - SILICONE MASK SIZE 4
PAL32070 - SILICONE MASK SIZE 5



DISPOSABLE RUBBER MASKS

Anatomic ventilation mask with blue rubber rim. Transparent shell allowing patient's monitoring during ventilation.

PAL32250 - RUBBER MASK SIZE 3
PAL32254 - RUBBER MASK SIZE 4
PAL32258 - RUBBER MASK SIZE 5



DISPOSABLE PVC MASKS

Disposable ventilation mask with turning standard connection (22 mm). Shell shaped transparent plastic material and border/cushion in soft inflated plastic.

PAL32150 - PVC MASK SIZE 0
PAL32154 - PVC MASK SIZE 1
PAL32158 - PVC MASK SIZE 2
PAL32162 - PVC MASK SIZE 3
PAL32166 - PVC MASK SIZE 4
PAL32170 - PVC MASK SIZE 5

FILTER AND RESERVOIR



PATIENT VALVE

Autoclavable patient valve for resuscitator bags OB 360° rotating connector. Adjustable safety over pressure valve.

PAL32370 - AUTOCLAVABLE PATIENT VALVE
PAL32375 - MEMBRANE FOR POS. 15 VALVE



DISPOSABLE RESERVOIR WITH EXTERNAL VALVE

Disposable reservoir for resuscitator bags OB with external valve.

PAL32350 - 2500 ml
PAL32356 - 900 ml



PATIENT VALVE WITH PAL34370 PEEP VALVE CONNECTION

Autoclavable patient valve for OB resuscitation bags 360° spinning connection. Adjustable over-pressure security valve. Equipped with PAL34371 Peep valve connection.

PAL34370 - PATIENT VALVE WITH PAL34370 PEEP VALVE CONNECTION



DISPOSABLE RESERVOIR WITH INTEGRATED VALVE

Disposable reservoir for resuscitator bags OB with integrated valve.

PAL34350 - 2500 ml
PAL34356 - 900 ml



PEEP VALVE - PAL34371

Peep Valve manufactured for patient valve PAL34370. Built in polycarbonate, autoclavable with steam for 30 cycles. Adjustable from 2,5 to 20 cmH₂O. It cannot be used on valves of other types or brands.



DISPOSABLE FILTER - 3505879

Disposable filter for resuscitator bags and automatic respirators. Sterile package.

VENTILATION SETS



REANIMATION BAG

Reanimation bag with shoulder strap for transport, featuring: resuscitator bag with turning connection security valve, 2 ventilation masks, reservoir, connecting tube for oxygen, 5 Guedel airways.

Dimensions: about 28x15x16 cm.

Weight: about 1 kg.

VAL40910 - REANIMATION BAG - SILICONE

VAL40914 - REANIMATION BAG - RUBBER

VAL40916 - REANIMATION BAG - PVC



RESUSCITATION CASE SILICONE OB - VAL40902

ABS resuscitation case with wall anchorage system. 1 autoclavable silicone resuscitator bag. Reservoir and oxygen tube. 3 silicone transparent ventilation masks. 5 Guedel airways. Dimensions: 39,5x13,5x26(h) cm. Weight: about 2,3 kg.

VAL40902 - SILICONE RESUSCITATION CASE

VAL40904 - RUBBER RESUSCITATION CASE

VAL40905 - PVC RESUSCITATION CASE



PEDIATRIC/NEWBORN RESUSCITATOR BAGS SETS

Set of 2 rubber resuscitator bags (newborn 240ml, pediatric 500ml) featuring: 3 ventilation masks, 1 reservoir, 1 oxygen tube, 3 Guedel airways.

Dimension: about 25x25x8(h) cm.

VAL44300 - RUBBER SET*



RELIABILITY

Diagnostic devices selected from the best manufacturers available on the market. All the devices we market are scrupulously checked according to documented procedures to guarantee quality and reliability.

STETHOSCOPES



SINGLE HEAD STETHOSCOPE

Compact and lightweight blue flat-head stethoscope. “Y” tube made of a single piece (PVC material). Very sensitive membrane with excellent tone transmission. CE marked according to the recent legal provisions concerning medical devices. Risk class according to Directive 93/42/CEE: I. The device complies with the REACH regulation in force.

DIA01012P - STETO. SINGLE HEAD PURPLE
DIA01012PK - STETO. SINGLE HEAD PINK
DIA01012Y - STETO. SINGLE HEAD YELLOW



SINGLE HEAD SATINED MASTER STET- DIA01020B

Professional single-head stethoscope. Highly sensitive level of auscultation. Head made of galvanized aluminum alloy with satin finish that guarantees lightness and durability. Excellent sound transmission thanks to the soft and flexible PVC tube. The package includes a pair of rigid cups, a soft pair and a membrane. Risk class according to Directive 93/42/CEE: I. The device complies with the REACH regulation in force.



SINGLE HEAD STETHOSCOPE - HS30A

Flat head stethoscope, compact and lightweight. Single binaural “Y” tube covered in molded PVC black colour. Sensitive membrane for a good sound transmission. Complies with Directive 93/42/EEC and subsequent amendments and supplements. Class I CE medical device.

SPHYGMOMANOMETERS



ANEROID SPHYGMOMANOMETER - DIA02003RE

Red aneroid sphygmomanometer. High visibility pressure gauge with graduated scale up to 300 mm/Hg. Nylon cuff with “D” ring and Velcro closing (allows a rapid and fast fastening). Separate inflation bulb with pressure regulation and release valve. Equipped with transport case. Risk class I m, according to MDD 93/42/CEE (Measuring instrument). The device complies with the REACH regulation in force.



PEDIATRIC PALMAR SPHYGMOMANOMETER - DIA100070

Aneroid palmar sphygmomanometer. Big manometer with graduated scale up to 300 mm/Hg. Complete with 2 cuffs, newborn and paediatric in fantasy fabric material with velcro closing. One-handed operations allowed, accurate regulation valve. Supplied with fantasy fabric case. Equipped with transport case. Risk class I m, according to MDD 93/42/CEE. The device complies with the REACH regulation in force.



SPHYGMOMANOMETER WITH STETHOSCOPE - HS-50C

Aneroid sphygmomanometer complete with integrated black stethoscope. Greater flexibility and ease of use (thanks to the integrated stethoscope). High visibility pressure gauge with graduated scale up to 300 mm/Hg. Bracelet with “D” ring, in black nylon with Velcro closure. Separate inflation pump with pressure relief and adjustment valve. Equipped with transport case. Measurement accuracy: ± 3 mmHg. Risk class according to Directive 93/42/CEE: I m. The device complies with the REACH regulation in force.

SPHYGMOMANOMETERS

LCD DIGITAL SPHYGMOMANOMETER - DM491



Automatic desktop digital sphygmomanometer. Innovative design with easy-to-read LCD display that simultaneously shows systolic and diastolic blood pressure, and heart rate. Display of date and time and detection of irregular heartbeat IHD (Arrhythmia). The device can store two sets of measurements for two patients (60 per patient). Oscillometric reading method according to international WHO standard. Automatic inflation cuff with Velcro closure. Pressure measuring range: 0 to 300 mmHg. Heart rate measurement range: 30 to 180 BPM. Accuracy: ± 3 % for pressure and ± 5 % for heartbeat. Power supply: 4 batteries (AA) 1.5 V, included. Autonomy: about 2 months with 3 measurements per day. Class IIa medical device. Dimensions: 148mm x 100mm x 56mm approx. Total weight: 382 gr (without batteries).

TABLE ANEROID SPHYGMOMANOMETER - DIA02009



Circular wall and/or table aneroid sphygmomanometer (diameter 170mm). Scale from 0 to 300 mm/Hg. Velcro-closing cuff and inflation bulb, spiral connecting tube. Wall anchor support with basket bracelet holder. Risk class I m, according to directive MDD 93/42/CEE. The device complies with the REACH regulation in force.

OMRON SPHYGMOMANOMETERS



OMRON M2 SPHYGMOMANOMETER - OMM2EC

Digital Automatic Blood Pressure Monitor. Compact, fully automatic, user-friendly and quick. Measures blood pressure and pulse rate. “IntelliSense” technology for comfortable controlled inflation. Detects irregular heartbeats. Memory for 30 sets. Digital LCD display. Power source: 4 AA batteries 1,5 V (included).



OMRON M3 SPHYGMOMANOMETER - OMM3IWIN

Digital arm blood pressure monitor. Compact, fully automatic, simple and fast. Comfortable controlled inflation with “Intelli Sense” technology. Precise measurement with “Intelli Wrap” comfort cuff. Indicator of correct armband positioning. Movement indicator. Detection of irregular heartbeats. Display of the average value of the last 3 measurements. Memory capacity of 60 measurements with date and time for each user (1 and 2). Power supply: 4 AA batteries (included). Class IIa medical device in compliance with MDD/93/42/EEC Directive.



OMRON M7 INTELLI IT SPHYGMOMANOMETER - OMM7ITIN

Digital arm blood pressure monitor. Compact, fully automatic, simple and fast. Automatic inflation through electric pump. Detection of irregular heartbeats. Memory capacity of 100 measurements per user. Digital LCD display. Power supply: 4 AA batteries (included) or mains adapter optionally. Class IIa medical device in compliance with MDD/93/42/EEC Directive as amended.

PULSE OXIMETERS



ADULT FINGER PULSE OXIMETER - FS20C-NS

Device for measuring blood oxygen saturation. Compact and lightweight, equipped with colour LED display. Detects oxygen saturation (SpO₂) and pulse rate. Battery level indicator. Without acoustic alarms. Provided with 2 AAA batteries. Dimensions: about 58x33x36mm. Maximum weight: 60g (with batteries). Class IIb CE medical device (93/42/CEE).

FS20C-NS - ADULT FINGER PULSE OXIMETER
FS-01 - CASE FOR FINGER OXIMETER



HANDHELD PULSE OXIMETER - HS10A

Lightweight for carrying and Easy-To-Use. Silicon rubber shell protection and stable metal bracket for table usage. Adult, pediatric and neonate oximeter probe available. CE medical device. Isolation class: II. Protection against electric shock: type BF – applied part. Protection against the penetration of fluids and substances: IPX1. Visual and acoustic alarm available. Big size 2.4" LCD display with SPO₂, pulse ratio, pulse bar-graph and plethysmogram. SPO₂ resolution 1 %. Low Battery voltage indicator. Refreshing rate: 1 sec. The unit can work with alkaline batteries (2000 mAh) or rechargeable Lion battery pack (1600 mAh).

HS10A - HANDHELD PULSE OXIMETER

DIAGNOSTIC PENLIGHTS



OB DIAGNOSTIC PENLIGHT

Rechargeable diagnostic bright penlight. Built in plastic material.
Pocket clip. Complete with 2 batteries AAA.

DIA101010 - OB DIAGNOSTIC PENLIGHT
DIA101020 - PEN LIGHT CLIP



DIAGNOSTIC METAL PENLIGHT

Rechargeable diagnostic bright penlight. Built in metal material.
Pocket clip. Complete with two batteries AAA.

DIA101011 - DIAGNOSTIC METAL PENLIGHT



DIAGNOSTIC PENLIGHT WITH LED

Rechargeable diagnostic bright penlight with LED light. Built in metal material. Pocket clip.
Complete with 2 batteries AAA.

DIA101012 - DIAGNOSTIC PENLIGHT WITH LED



SET 6 DIAGNOSTIC PENLIGHTS

Set of 6 disposable diagnostic pens with filament lamps. Made of plastic material and water resistant. To turn on the light simply press the clip on the pen itself. Disposable device, battery not replaceable.

DIA05004 - SET 6 PCS DISPOSABLE
DIAGNOSTIC PENLIGHTS

THERMOMETERS



INFRARED DIGITAL THERMOMETER - FT100B

LCD display with 3½ figures. Temperature unit of measurement: Celsius degrees (°C) and Fahrenheit (°F). Measurement range: 32.0 – 42.9 °C. Saves last 10 performed measurements. Powered by not rechargeable battery, PC 2032 type (1pc). The batteries last about 2500 uses. Size: 105x35x35mm. Net weight: 36gr. Acoustic alarm: beep sound for ON-OFF. Complete with protective plastic case.



DIGITAL THERMOMETER - RM300

Measurement of body temperature by oral, axillary or rectal contact. Measurement in 60 seconds. Units of measurement in degrees Celsius (°C) and Fahrenheit (°F). Sensed temperature range: 32°C to 42.9°C. LCD display with 3 displayed digits. Water resistant tip (IP27 degree of protection). 1.5 V alkaline battery. Autonomy of approx. 200 hours. Keeps track of the last measurement. Class IIa EC medical device in compliance with MDD/93/42/EEC and subsequent amendments and additions.



MERCURY FREE THERMOMETER- TER174

Clinical mercury-free thermometer for body temperature detection. The liquid used is a "Gallio/Indio" alloy combination, totally non-toxic and safe for the environment. CE medical device. Detection range: 35-42° C, precision: 0,10° C. Dimensions: about 125x12x8 mm.

DIAGNOSTICS DEVICES AND ACCESSORIES



PRESSURE INFUSER WITH OB MANOMETER - DIA 100390

Pressure infuser for medicaments contained in plastic material sachets. Complete with manometer for infusion monitoring and mono-ball with deflating valve for accurate dosage. Suitable for 500 ml infusion. CE medical device.



BLACK OPHTHALMOSCOPE - DIA 100415

Fiber optic 2,5V bulb. Ophthalmoscope: head with screw blocking; corrective lenses from 0 to +20 and 0 to -20 dioptres with 5 different diaphragm apertures. Plastic/metal handle with batteries; working with 2 batteries type AA (not included). Complete with shockproof plastic case. CE medical device.



WOODEN TONGUE DEPRESSOR - 90554

Wooden tongue depressor. Separately packaged. Dimensions: 150x17 mm.



MINI OTOSCOPE WITH CONVENTIONAL LIGHT - DMV540B

High quality blue otoscope for pocket and doctor's bag. Unobstructed view, unit equipped with halogen bulb. Robust and light weight plastic with metal inserts and 3x rotating magnifying glass. Handle requires 2 batteries (AA) and it is equipped with grooves for a better grip. Speculum input is rapid and stable. Provided with bag and with 4 additional reusable specula.

DIAGNOSTICS DEVICES AND ACCESSORIES



ECG ELECTRODES

E.C.G. electrodes. Base made of foam material. Conductive paste. Optimal detection of the tracing.

Adult: diameter 55 mm.

Pediatric: diameter 43 mm.

50 pcs package.

F9060 - ADULT ECG ELECTRODES 50 pcs


F9060P - PEDIATRIC ECG ELECTRODES 50 pcs



ADULT PADS - F7956W

PADs for adults.

Suitable for defibrillator OB Life Support 300. PADs with blue connector should be used for patients, whose weight is more than 25 kg.

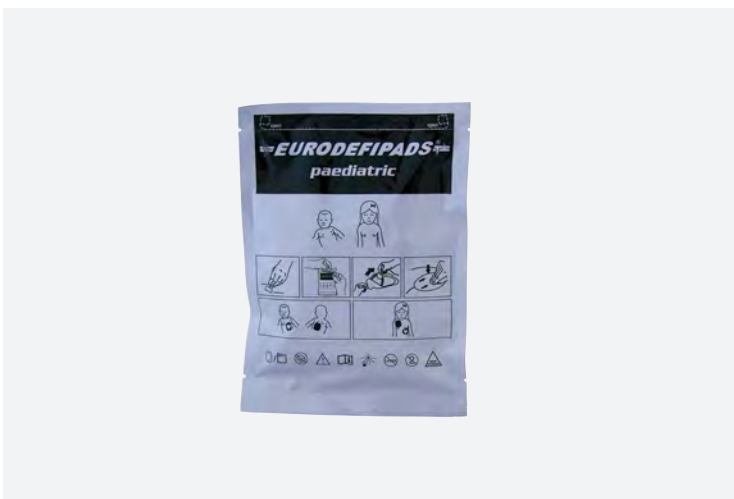

PEDIATRIC PADS - F7956P

PADs for children.

Suitable for defibrillator OB Life Support 300. PADs with yellow connector should be used for children, whose weight is less than 25 kg (younger than 8 years of age).



REANIMATION

We have a range of products for the most critical situations: medical devices for reanimation and specific kits. Our goal of guaranteeing reliability and safety for patients and users is achieved through constant quality controls.

CONVENTIONAL LIGHT LARYNGOSCOPE



MC-INTOSH LARYNGOSCOPE

Mc Intosh blades for conventional light laryngoscope, made in steel. Different sizes available. Sets with autoclavable laryngoscope available with 3 blades (sizes from 1 to 3) or 4 blades (sizes from 0 to 3). Sets are supplied with a practical roll case.

RIA53000 - 3 BLADES SET MC-INTOSH CONVENTIONAL LIGHT
RIA53004 - 4 BLADES SET MC-INTOSH CONVENTIONAL LIGHT
RIA53020 - MC-INTOSH CONVENTIONAL LIGHT BLADE SIZE 0
RIA53022 - MC-INTOSH CONVENTIONAL LIGHT BLADE SIZE 1
RIA53024 - MC-INTOSH CONVENTIONAL LIGHT BLADE SIZE 2
RIA53026 - MC-INTOSH CONVENTIONAL LIGHT BLADE SIZE 3
RIA53028 - MC-INTOSH CONVENTIONAL LIGHT BLADE SIZE 4



MILLER LARYNGOSCOPE

Conventional light metal autoclavable laryngoscope with straight Miller blades in steel for conventional light laryngoscope. Blades are available in different sizes. Available in set with 3 blades (sizes from 0 to 2) and a practical roll case.

RIA53100 - 3 BLADES SET MILLER CONVENTIONAL LIGHT
RIA53120 - MILLER CONVENTIONAL LIGHT BLADE SIZE 0
RIA53122 - MILLER CONVENTIONAL LIGHT BLADE SIZE 1
RIA53124 - MILLER CONVENTIONAL LIGHT BLADE SIZE 2
RIA53126 - MILLER CONVENTIONAL LIGHT BLADE SIZE 3
RIA53128 - MILLER CONVENTIONAL LIGHT BLADE SIZE 4

FIBER OPTIC LARYNGOSCOPE



MC-INTOSH LARYNGOSCOPE

Mc Intosh fiber optic blades made in steel for fiber optic laryngoscope. Different sizes available. 2,5V halogen lamp set inside the blade for an easy sterilization. Sets with laryngoscope available with 3 blades (sizes from 1 to 3) or 4 blades (sizes from 0 to 3). Sets are supplied with a practical roll case.

RIA53200 - 3 BLADES SET MC-INTOSH FIBER OPTIC
RIA53204 - 4 BLADES SET MC-INTOSH FIBER OPTIC
RIA53220 - BLADE MC-INTOSH FIBER OPTIC SIZE 0
RIA53222 - BLADE MC-INTOSH FIBER OPTIC SIZE 1
RIA53224 - BLADE MC-INTOSH FIBER OPTIC SIZE 2
RIA53226 - BLADE MC-INTOSH FIBER OPTIC SIZE 3
RIA53228 - BLADE MC-INTOSH FIBER OPTIC SIZE 4
RIA53230 - BLADE MC-INTOSH FIBER OPTIC SIZE 5



MILLER LARYNGOSCOPE

Miller fiber optic blades made in steel for fiber optic laryngoscope. Different sizes available. 2,5V halogen lamp set inside the blade for an easy sterilization. Sets with laryngoscope available with 2 blades (sizes from 0 to 1) or 3 blades (sizes from 0 to 2). Sets are supplied with a practical roll case.

RIA53300 - 2 BLADES SET MILLER FIBER OPTIC
RIA53304 - 3 BLADES SET MILLER FIBER OPTIC
RIA53320 - BLADE MILLER FIBER OPTIC SIZE 0
RIA53322 - BLADE MILLER FIBER OPTIC SIZE 1
RIA53324 - BLADE MILLER FIBER OPTIC SIZE 2
RIA53326 - BLADE MILLER FIBER OPTIC SIZE 3

FIBER OPTIC LARYNGOSCOPE WITH DISPOSABLE BLADES



MC-INTOSH DISPOSABLE BLADES

Disposable fiber optic Mc Intosh laryngoscope blade, made in plastic material. Different sizes available.

Size 1: 92x15(h)

Size 2: 112x15(h)

Size 3: 132x25(h)

Size 4: 165x30(h)

RIA53444 - DISPOSABLE MC INTOSH BLADE SIZE 1

RIA53446 - DISPOSABLE MC INTOSH BLADE SIZE 2

RIA53448 - DISPOSABLE MC INTOSH BLADE SIZE 3

RIA53450 - DISPOSABLE MC INTOSH BLADE SIZE 4



MILLER DISPOSABLE BLADES

Disposable fiber optic Miller laryngoscope blade, made in plastic material. Different sizes available.

Size 0: 75x8(h)

Size 1: 100x8(h)

Size 2: 132x9(h)

Size 3: 172x12(h)

RIA53432 - DISPOSABLE MILLER BLADE SIZE 0

RIA53434 - DISPOSABLE MILLER BLADE SIZE 1

RIA53436 - DISPOSABLE MILLER BLADE SIZE 2

RIA53438 - DISPOSABLE MILLER BLADE SIZE 3



3 FIBER OPTIC DISPOSABLE BLADES SET - RIA53410

Fiber optic laryngoscope complete with 3 Mc Intosh blades (2-3-4). Made in plastic material. Autoclavable handle. Supplied with a practical roll case.



6 FIBER OPTIC DISPOSABLE BLADES SET - RIA53400

Fiber optic laryngoscope complete with 6 blades (3 Mc Intosh, 3 Miller) and handle. Supplied with a practical case. Dimensions: about 28x4,5x25 cm.

LARYNGOSCOPE ACCESSORIES



CONVENT. LIGHT LARYNG. HANDLE

Conventional light laryngoscope handle made in stainless steel.
Dimensions: Ø 30x148 mm.

RIA53030 - CL L. HANDLE
RIA53182 - CL L. BLADE BULB



FIBER OPTIC LARYNG. HANDLE

Fiber optic laryngoscope handle made in stainless steel.
Dimensions: Ø 30x152 mm.

RIA53232 - FO L. HANDLE
RIA53380 - FO L. BLADE BULB



SMALL CONVENT. LIGHT LAR. HANDLE

Small conventional light laryngoscope handle made in stainless steel.
Dimensions: Ø 20x155 mm.

RIA53130 - SMALL CL L. HANDLE
RIA53180 - CL L. SMALL BLADE BULB



SMALL FIBER OPTIC LARYNG. HANDLE

Small fiber optic laryngoscope handle made in stainless steel.
Dimensions: Ø 20x152 mm.

RIA53330 - SMALL FO L. HANDLE
RIA53382 - FO L. SMALL BLADE BULB



DISPOSABLE FIBER OPTIC L. HANDLE

Fiber optic disposable laryngoscope handle, autoclavable.
Dimensions: Ø 30x154 mm.

RIA53454 - DISPOSABLE FO L. HANDLE
RIA53458 - FO DISPOSABLE L. BLADE BULB



ROLL-CASE FOR LARYNGOSCHOPE

Roll-case for laryngoscope - red. 5 inner compartments. Made in washable cloth with inner stuffing.
Dimensions : about 19x34 cm.

BOR84041 - ROLL CASE FOR LARYNGOSCHOPE

AIRWAYS AND CATHETER MOUNT



GUEDEL AIRWAY NO. 000

Pink - 40 mm
Ref.: 127141001004



GUEDEL AIRWAY NO. 00

Blue - 50 mm
Ref.: 127141001006



GUEDEL AIRWAY NO. 0

Black - 60 mm
Ref.: 127141001008



GUEDEL AIRWAY NO. 1

White - 70 mm
Ref.: 127141001010



GUEDEL AIRWAY NO. 2

Green - 80 mm
Ref.: 127141001012



GUEDEL AIRWAY NO. 3

Yellow - 90 mm
Ref.: 127141001014



GUEDEL AIRWAY NO. 4

Red - 100 mm
Ref.: 127141001016



GUEDEL AIRWAY NO. 5

Light blue - 110 mm
Ref.: 127141001018

RIA50290 - SET OF 5 GUEDEL AIRWAYS
(0/1/2/3/4)



WENDL AIRWAY

Wendl nasopharyngeal airway, with ending pipe ring. Supplied singularly in sterile package. Different sizes available.

185200CH26 - WENDL AIRWAY CH 26-8,7 mm
185200CH28 - WENDL AIRWAY CH 28-9,3 mm
185200CH30 - WENDL AIRWAY CH 30-10 mm



CATHETER MOUNT - 1515TG

Catheter Mount, lightweight, transparent, disposable and at low compliance. Supplied with perfect capacity connectors. Extremely flexible. Allows up to 360° angulation without any modification to inside volume, crushproof. Connectors in compliance to ISO norms. Dimensions: 20 cm.

REANIMATION DEVICES



STERILE ENDOTRACHEAL TUBES

Cuffed endotracheal tubes in plastic material. Unidirectional valve, standard connector in accordance to ISO norms. Reference notches, x-ray opaque, non-traumatic tip. Supplied in sterile package. Different sizes available.

11814102003	- CH12 3,0MM
11814102004	- CH14 3,5MM
11814102005	- CH16 4,0MM
11814102006	- CH18 4,5MM
11814102007	- CH20 5,0MM
11814102008	- CH22 5,5MM
11814102009	- CH24 6,0MM
11814102010	- CH26 6,5MM
11814102011	- CH28 7,0MM
11814102012	- CH30 7,5MM
11814102013	- CH32 8,0MM
11814102014	- CH34 8,5MM



THOMAS HOLDER

Safety endotracheal tubes blocking system. Screw system, integrating a no-lockjaw surface which protects teeth and gums and avoids endotracheal tube crushing. It avoids use of plaster.

05-1650	- ADULT THOMAS HOLDER
05-1651	- PEDIATRIC THOMAS HOLDER

REANIMATION DEVICES



LIFE KEY MASK KEY-RING

Key-ring with mouth to mouth ventilation mask. Elastics for ear fixing. One way valve for a safe ventilation. Singularly packaged . Non-sterile.

RIA50248 - BLUE LIFE KEY MASK KEY-RING
FS-102 - LIFE KEY MASK



QUICKTRACH

Urgent coniotomy sterile set quicktrach type. Complete with 1 plastic canula with intubation internal-blade stylet with diameter of 4 mm for adults and 2 mm for child. Complete with ISO joint for the connection to resuscitation bag. Quicktrach features a sharp conical-tip intubation stylet that avoids windpipe incision using classic scalpel and a "lock" that avoids to penetrate too much in depth. 2 lateral wings allowing airways safe fixing.

VBM30-04-0021 - CHILD QUICKTRACH
VBM30-04-0041 - ADULT QUICKTRACH



REMASK VENTILATION MASK - RIA50104

Remask transparent ventilation Mask, with pneumatic border. Mask adjustable to all kind of resuscitator bag. Disposable one-way safety valve. Blocking elastic. Connecting juncture to oxygen source. Standard connector 22mm. Supplied in a practical yellow case. CE certified in class IIa.



INTUBATION STYLETS

Flexible and coated intubation stylets. Different sizes available.

RIA52213 - INTUBATION STYLET FR6 - 2,0 MM
RIA52214 - INTUBATION STYLET FR10 - 3,3 MM
RIA52215 - INTUBATION STYLET FR14 - 4,7 MM

INTUBATION BAG - CONVENTIONAL LIGHT



INTUBATION BAG COMPLETE - RIA50820

Intubation bag complete of: 1 conventional light laryngoscope with 3 blades (1-2-3), 1 Magill forceps, 1 Pean forceps, 1 scissors, 5 endotracheal tubes, 3 Guedel airways, 1 syringe, 1 mouth opener, 1 plaster roll, 2 Wendl airways, 1 stylet.

Dimensions: about 32x22x7,5(h) cm.

Weight: 1,2 kg.



PEDIATRIC INTUBATION BAG COMPLETE - RIA50816

Intubation bag complete of: 1 pediatric conventional light laryngoscope with 3 blades (0-1-2), 1 Magill forceps, 1 Pean forceps, 1 scissors, 4 endotracheal tubes, 3 Guedel airways, 1 syringe, 1 mouth opener, 1 plaster roll, 2 Wendl airways, 1 stylet.

Dimensions: about 32x22x7,5(h) cm.

Weight: 1 kg.



INTUBATION SET IN POUCH - VAL44410

Complete intubation set in pouch with transparent frontal, featuring: 1 laryngoscope pediatric handle, 1 C.L. Mc Intosh blade size 0, 1 C.L. Miller blade size 1, 1 C.L. Miller blade size 0, 1 Magill forceps, size 16, 1 Pean forceps, 4 endotracheal tubes, 4 Guedel airways sizes 000/00/0/1, 1 intubation stylet 1 plaster roll.

Dimensions: about 27x12x9(h) cm.

Weight: 600 g.



PEDIATRIC INTUBATION SET IN POUCH - VAL44390

Complete intubation set in pouch with transparent frontal, featuring: 1 laryngoscope with 3 conventional light blades, 1 Magill forceps, 1 Pean forceps, 1 scissors, 5 endotracheal tubes, 3 Guedel airways, 1 syringe, 1 mouth opener, 1 plaster roll, 2 Wendl airways, 1 Intubation stylet.

Dimensions: about 27x12x9h cm.

Weight: 1kg.

INTUBATION BAG - FIBER OPTIC



INTUBATION BAG COMPLETE - RIA50822

Intubation bag complete of: 1 fiber optic laryngoscope with 3 blades (1-2-3), 1 Magill forceps, 1 Pean forceps, 1 scissors, 5 endotracheal tubes, 3 Guedel airways, 1 syringe, 1 mouth opener, 1 plaster roll, 2 Wendl airways, 1 stylet.

Dimensions: about 32x22x7,5(h) cm.

Weight: 1,3 kg.



PEDIATRIC INTUBATION BAG COMPLETE - RIA50825

Intubation bag complete of: 1 fiber optic laryngoscope with 3 blades (0-1-2), 1 Magill forceps, 1 Pean forceps, 1 scissors, 4 endotracheal tubes, 3 Guedel airways, 1 syringe, 1 mouth opener, 1 plaster roll, 2 Wendl airways, 1 stylet.

Dimensions: about 32x22x7,5(h) cm.

Weight: 1 kg.



DISPOSABLE INTUBATION BAG COMPLETE- RIA50826

Intubation bag complete of: 1 fiber optic laryngoscope with 6 blades, 1 Magill forceps, 1 Pean forceps, 1 scissors, 5 endotracheal tubes, 3 Guedel airways, 1 syringe, 1 mouth opener, 1 plaster roll, 2 Wendl airways, 1 stylet.

Dimensions: about 32x22x7,5(h) cm.

Weight: 1kg.



SURGICAL MEDICATION AND INSTRUMENTS

We distribute surgical medication and instruments that are subjected to quality controls before being placed on the market.

SPECIFIC MEDICATION KITS

We create specific kits to deal with emergencies of various entities and nature, for domestic, non-professional and professional use.

FIRST AID CASES

We provide first aid boxes that meet the requirements of workplace safety.

GAUZE COMPRESSES



STERILE GAUZE COMPRESSES 2 PCS

Sterile medication gauze compresses in an easy-to open package, two 8-ply 10x10 compresses, pure cotton tit.12/8. Indicated for emergency and first aid activities. 2 pcs package. Ref.: MED96001.



NON WOVEN STERILE GAUZE COMPRESSES 2 PCS

Sterile non woven compresses. Absorbing, non-adherent and resistant, two 8-ply. 10x10 compresses. 2 pcs package. Ref.: MED96101.



STERILE GAUZE COMPRESS

Sterile medication gauze compresses in an easy-to open package, 1 18x40 cm compress, pure cotton tit.12/8. Indicated for emergency and first aid activities. Ref.: MED96008.



NON WOVEN SELF-ADHESIVE

Non woven self-adhesive roll for medication fixing on all body surfaces. Dimensions: 10 m x15 cm. Ref.: 35CER334.



1 KG NON-STERILE GAUZE COMPRESS

Non-sterile medication gauze, pure cotton. 1 kg package. Different sizes available.

61610 - 1 kg GAUZE 10x10 cm
61645 - 1 kg GAUZE 15x15 cm
61670 - 1 kg GAUZE 20x20 cm
61710 - 1 kg GAUZE 30x30 cm




NON-STERILE GAUZE COMPRESSES – 100PCS

12-layer, non-sterile gauze compresses made of pure cotton. Package: 100 pcs. Different sizes available.

61050 - 100 pcs 5x5 cm
61051 - 100 pcs 7,5x7,5 cm
61052 - 100 pcs 10x10 cm

ELASTIC TUBULAR NET, BANDAGES AND COTTON

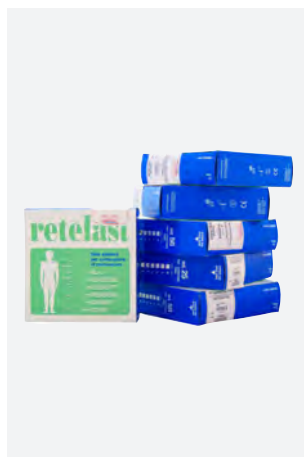

ELASTIC TUBULAR NET

Elastic tubular net for medications on every body's part. Easy and quick application.

75412 - SIZE 1 - 50 m
75414 - SIZE 2 - 50 m
75416 - SIZE 3 - 50 m
75418 - SIZE 4 - 50 m
75420 - SIZE 5 - 50 m
75422 - SIZE 6 - 25 m



ELASTIC TUBULAR NET

Elastic tubular net for medications. Easy and quick application. Length: about 40cm. The size is indicated for hand, wrist, foot, ankle.

75466 - ELASTIC TUBULAR NET MEDIUM



ELASTIC MEDICATION BANDAGE

Elastic medication bandage. Singularly packaged. Dimensions: about 4 m x 10 cm

MED96280 - ELASTIC BANDAGE



EDGED BANDAGE

OB flexible medication bandage. Singularly packaged in sealed wrapping. Different sizes available.

FS-3570 - BANDAGE 5 cm
FS-3572 - BANDAGE 10 cm
FS-3554 - BANDAGE 8 cm



COHESIVE BANDAGE

Self-fixing elastic bandage, lightweight and soft. Self adherent thanks to its atural rubber latex micro particles. Does not adhere to the skin and it is highly transpiring. Different sizes available, lenght 4 m.

370 - COHESIVE BAND. 4 cm
371 - COHESIVE BAND. 6 cm
372 - COHESIVE BAND. 8 cm
373 - COHESIVE BAND. 10 cm



COTTON

Absorbent cotton wool 20 or 50 g package. Not a medical device.

COT104 - COTTON 50 g
COT279 - COTTON 20 g

SOFT NEXT COHESIVE BANDAGE AND DETECTABLE DRESSINGS

SOFT NEXT COHESIVE BANDAGE

Ideal for covering wounds especially at critical points, such as hands fingers. It is used directly on the wound. It is mouldable, autostick, elastic, water resistant and latex free. It absorbs the blood and does not stick to the skin, the hair and the wound itself. Skin or blue colour, different sizes available.

CER050 - SOFT NEXT 3x450 cm SKIN
CER051 - SOFT NEXT 3x450 cm BLUE
CER059 - SOFT NEXT 6x100 cm SKIN
CER064 - SOFT NEXT 6x100 cm BLUE
CER063 - SOFT NEXT 6x450 cm SKIN
CER061 - SOFT NEXT 6x450 cm BLUE



DETECTABLE BLUE DRESSINGS - CER049

Pack of 100 detectable sterile blue dressings, 5 assorted sizes. The metal band fitted in the internal part of the dressing is easy to find by the metal detectors. Contains: 20 pcs 9x72 mm, 40 pcs 30x72 mm, 32 pcs 55x18 mm, 4 pcs butterfly, 4 pcs 72x50 mm.



EXTRA LONG 50PCS BLUE DETECTABLE DRESSING - CER143

The dressing length helps you to wrap your finger more than once, protecting the eventual cut or graze and avoiding accidental detachments, as it may happen for the smaller dressings. Dimensions: 18x2 cm. Pack of 50pcs.



TAPE, PLASTER AND SKIN SAVER

NEURO-MUSCULAR KINESIOLOGY TAPE

Neuro-muscular tape, Thera-Band® Kinesiology Tape eliminates the risks of incorrect application thanks to the exclusive XactStretch™ Technology: our guide indicators allow you to always apply the tape with the correct tension. Relieves pain. Provides support for muscles and joints. Assists the treatment process of the edema. 98% Cotton, 2% Sapdex. Latex Free. Dimensions: 5 cm x 5 m.

305464 - KINESIOLOGY TAPE BLACK/GREY
305465 - KINESIOLOGY TAPE WHITE/BLACK
305467 - KINESIOLOGY TAPE BLUE
305468 - KINESIOLOGY TAPE PINK/WHITE
305469 - KINESIOLOGY TAPE RED/BLACK
305470 - KINESIOLOGY TAPE GREEN/YELLOW



TAPE PLASTER

Sport tape in white viscose with zinc oxide adhesive, suitable for rigid bandaging. Grainy fray-proof borders for an easy snatch. Available in two dimensions for a length of 14 m.

35CER692 - TAPE 3,8 cm
35CER694 - TAPE 5 cm



ELASTICATED STRIP PLASTER

Medicated strip plaster in elasticated cotton material, with transpiring gauze pad. Available in two dimensions for a width of 6 cm.

70051 - STRIP PLASTER 1 m
70054 - STRIP PLASTER 5 m



SKIN SAVER - MED97049

Thin and light layer made of polyurethane foam. Apply before bandaging. Dimensions: cm 7x27,5 m.



POLYETHYLENE STRIP PLASTER - 70104

Skin coloured polyethylene medicated strip plaster. Waterproof, elastic, transpiring and hypoallergenic. High absorbing viscous gauze pad, does not stick to the wound. Dimensions: 5 m x 6 cm. Singly packaged.



PLASTERS AND PLASTER ROLL

NON WOVEN PLASTER ROLL



Non woven plaster roll.
Available in two different
sizes with a lenght of 9,14
m.

00181 - ROLL 2,5 cm
00183 - ROLL 5 cm

ASSORTED PLASTERS - 400187



Breathable sterile classic
plaster assorted in 4
different sizes. Class
Is medical device in
compliance with Directive
93/42/EEC. Package
of 10 pieces. Singularly
packaged.


SILK PLASTER ROLL



Adhesive silk plaster roll
Available in two different
sizes with a lenght of 9,14
m.

35CER196 - SILK ROLL 2,5 cm
02023173100000 - SOFFIX
SILK ROLL 5 cm

ASSORTED PLASTERS

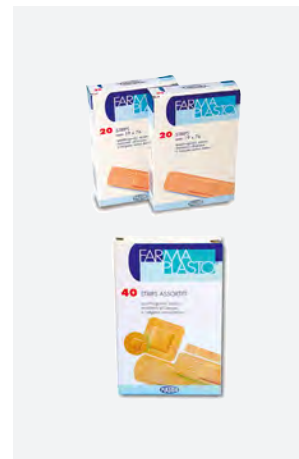

Medicated plasters with
gauze pad on adhesive
support. Singularly
packaged.

35CER212 - 20 PIECES
35CER216 - 40 PIECES

CLOTH PLASTER ROLL



Adhesive cloth plaster roll.
Available in two different
sizes with a lenght of 5 m.

02023132000000 -
CLOTH PLASTER 2,5 cm
02023133000000 -
CLOTH PLASTER 5 cm

CLOTH PLASTER ROLL - 400014



Adhesive cloth roll. Class
I non-sterile medical
device in compliance with
Directive 93/42/EEC. Single
pack.

PLASTERS



NON WOVEN PLASTER

Adhesive, sterile medication with central absorbing pad non adhering to wound
Hypoallergenic transpiring adhesive. Available in different sizes with a lenght of 10 cm.

35CER342 - PLASTER 8 cm
35CER346 - PLASTER 15 cm
35CER348 - PLASTER 20 cm



ADHESIVE STERILE EYE BANDAGE - 35CER314

Sterile medication ocular compress with a high absorbent power on adhesive support.
Hypoallergenic and transpiring.
For secure and comfortable application on the eye.



OMNISTRIP SUTURE PLASTER - 540684

Sterile suture plaster skin colour, coated with an hypoallergenic adhesive layer. Suitable for small wound non-traumatic closing. 10 sterile strips.



PLASTER FOR LINE-IN CANNULA NEEDLE - 35CER400

Non-woven plaster used to fasten and secure cannula needles. One package, 50 pcs. Size: 9x6 cm.

IMMOBILIZATION AND MEDICATION SHEET



TRIANGULAR IMMOBILIZATION SHEET

Triangular immobilization sheet, singularly packaged.

Available in cloth: 90x90x129 cm.

Available in non-woven: 90x90x127 cm.

FS-140 - CLOTH MATERIAL

FS-143 - NON-WOVEN MATERIAL



STERILE MEDICATION CLOTH

Sterile medication cloth, ready to use for physiological treatment of wounds and burns. Hypoallergenic, transpiring with high. Absorbing capability. Allows transpiration of wound or burn and absorbs without sticking. Provided in a single sterile wrap.

62050 - STERIL CLOTH 40x60 cm

62051 - STERIL CLOTH 60x80 cm

62052 - STERIL CLOTH 80x120 cm



HEMOSTATIC BANDAGE PAD - 1002001

Compressive bandaging for wounds and haemorrhages. It assures wound hemostasis thanks to the compression performed by elastic bandage on the hemostatic pad. The pad assures the hemostasis and moreover does not allow blood leakage. Supplied in sterile package. Dimensions: about 300x10 cm.

MEDICATION PACKAGES



MEDICAL PACKAGE D.I.N.

Sterile dressing package according to norm DIN 13151. Ready to use bandage and gauze. Hypoallergenic and with great transpiration power. Highly absorbent. Sterile packaging, immediately ready to use. Bandage length: 4 m. Different sizes available.

62000 - DIN MEDICATION PACKAGE S - 6 cm
62001 - DIN MEDICATION PACKAGE M - 8 cm
62002 - DIN MEDICATION PACKAGE L - 10 cm



ALUMINIUM MEDICATION PACKAGES, ALUMINIUM COMPRESS AND ALUMINIUM CLOTHS

Sterile dressing for wounds and burns: aluminized medication package. Hypoallergenic, transpiring, anti-bacteria, highly absorbent and incohesive. Composed of a layer of aluminum and a layer of gauze with high absorbency that greatly reduces the infection. The aluminum film, in direct contact with the wound or with the burn, allows the transudate or blood to pass through with the aluminum layer not adhering to the wound. Promotes epithelialization and scar tissue formation. Thanks to the great absorbent power, but not adherent, it facilitates the change of the dressing.

63050 - ALUMINIUM MEDICATION PACKAGE S
63051 - ALUMINIUM MEDICATION PACKAGE M
63052 - ALUMINIUM MEDICATION PACKAGE G
63003 - ALUMINIUM COMPRESS 10x10 cm
63100 - ALUMINIUM CLOTH 40x60cm
63101 - ALUMINIUM CLOTH 60x80 cm

DISPOSABLE GLOVES IN LATEX



STERILE LATEX

Chirurgical sterile gloves made of pure latex from natural rubber, with powder with anatomic shape. Enlarged backside, lubricated with corn starch and provided with ring beam. Singularly packaged with double wrapping made of medical paper 1 pair.

CLASSIC65 - SIZE 6,5
CLASSIC70 - SIZE 7,0
CLASSIC75 - SIZE 7,5
CLASSIC85 - SIZE 8,5
40400DM - GLOVES BOXES WALL BRACKET



STERILE LATEX SIZE L

1 pair of disposable sterile latex gloves.
Singularly packaged. Size: L.

01BML0S0G - STERILE
LATEX GLOVES SIZE L

TOURNIQUET



LATEX FREE

Quick release latex free tourniquet, one hand operated. Easy to use, ensuring an excellent hemostasis without skin and tissues laceration. Quick opening and progressive/adjustable haemostasis reduction.

MED98210 - RED
MED98212 - BLUE



FLAT TOURNIQUET 3 PCS - 500057

Set of 3 pcs flat nitrile tourniquets



LATEX FREE BABY - MED98208

Quick release latex free tourniquet, one hand operated. Easy to use, ensuring an excellent hemostasis without skin and tissues laceration. Quick opening and progressive/adjustable haemostasis reduction. Fantasy colour.



OB HEMOSTATIC BAND - MED98230

Tourniquet OB type Martin Wide surface mm 75 x 1000 mm for an excellent haemostasis. Suitable for huge wounds and amputations. Supplied with regulation non-traumatic buttons. Latex free.



FLAT TOURNIQUET - LAC190

Flat tourniquet. Singularly packaged. Latex free.



CAT COMBAT APPLICATION TOURNIQUET - CAT-007

Original Cat Combat Application Tourniquet haemostatic tourniquet, in nylon material with micro metrics regulation.

NEEDLES BOX



MULTI-SAFE MINI 300

Container suitable for needles and sharps. With safety locking system. Made of polypropylene. Volume: 0,3 L. Dimensions: 7,8x4,5x16 cm.

MED98609 - MULTI-SAFE MINI 300



MULTI-SAFE SANI 200

Container suitable for needles and sharps. With safety locking system. Made of polypropylene. Volume: 0,2 L. Dimensions: 7,3x4,5x13 cm.

MED98613 - MULTI-SAFE SANI 200



MULTI-SAFE VARIO 2000

Container suitable for needles and sharps. With safety locking system. Easy handling during lid assembly due to snap-on hooks. Clear fill level marking thanks to new label design. Compatible with standard trays. Made of polypropylene. Volume: 2,0 l. Dimensions: 95 x 240 x140 mm

MED98625 - MULTI-SAFE VARIO 2000



RIGIBOX EUROMATIC

Waste and needle box, specific for ambulances. With safety closure. Made of eco-friendly polypropylene with a red cap. Volume: 7 l. Size: about 17x17x 39 cm.

MED98631 - RIGIBOX EUROMATIC
MED98623 - WALL HOLDER FOR EUROMATIC



KONTAMED

Trash and needles box suitable for ambulance. Wall bracket for wall fixing. Complete of safety closure. Manufactured in yellow high resistant plastic material. Dimensions: about 170x75x225h mm.

2009067 - KONTAMED



KONTAMED MINI

Small trash and needles box suitable for emergency bags and cases. Complete of safety closure. Manufactured in yellow high resistant plastic material. Dimensions: 45x115h mm.

2009068 - KONTAMED MINI

SYRINGE AND NEEDLES

CANNULA NEEDLE



2 FEB ways intravenous catheter needle with pliable wings injection valve. Best convenient use, utmost fixing to the skin. Sterile, disposable, non-toxic, x-ray opacity, nonpyrogenic.

50AC2S14 - ORAN. 14-45 mm
50AC2S16 - GREY 16-45 mm
50AC2S18 - GREEN 18-45 mm
50AC2S20 - PINK 20-32 mm
50AC2S22 - L. BLUE 20-25 mm
50AC2S24 - YELL. 24-23 mm
35CER400 - PLASTER FOR
LINE-IN NEEDLE 50 PCS

BUTTERFLY NEEDLES



Sterile epicranial butterfly needle, for intravenous infusion. Indicated for patients with capillary fragility, for children, and in the case a thin atraumatic needle is needed. Easily flexible.

90185 - GREEN 21G
90187 - LIGHT BLUE 23 G

CUP CANNULA NEEDLE - 90227



Disposable sterile protective cap for needle canula.

3-WAY TAP - 90232



3-way tap with spinning luer lock.

STERILE SYRINGE



Sterile disposable syringes, singularly packaged.
1 Package, 100 pcs.

90152 - SYRINGE 2,5 cc
90153 - SYRINGE 5 cc
90154 - SYRINGE 10 cc
90155 - SYRINGE 20 cc

FLUIDS ADMINISTRATION SET - 01D193E2N



Infusion set for administration of solutions. (Not for blood or blood products). Sterile infusion set for delivery of liquids and fluids in intravenous therapy and medical application. With an antibacterial drip chambers, lines and tubing with roller clamp. Disposable with protective end caps and filter.

PROTECTIVE MASKS



ADULT SURGICAL MASK

Three-layer type IIR medical surgical mask with elastics to protect the nose and mouth. Disposable not sterile medical device. Non-woven fabric, without seams, with elastic bands applied without the use of adhesives. CE marking according to European Directive 93/42 CEE (D.L. 46/97). Risk class: I. Certified according to ISO 14683:2019 for medical facial masks.

FM003 - 10 PIECES

FPI002R - 50 PIECES



PEDIATRIC SURGICAL MASK - YX005

Surgical, non-sterile, disposable facial mask for paediatric use. Blue, with soft and light adjustable nose pad. Composed of 3 layers of TNT (non-woven fabric, latex free). Bacterial filtration efficiency (BFE) \geq 95%. Fluid resistant, dustproof, windproof, good breathability. Class I medical device compliant with Regulation (EU) 2017/745. Certified according to EN 14683:2019 for Type I medical use face masks. Bio-compatible in accordance with ISO 10993-1:2018. Pack of 10 pcs.



PROTECTIVE MASK FFP2

FFP2 (KN95) protective mask – not sterile, disposable intended for an adult user. Equipped with ear-loop and adjustable nose pad for light pressure on the skin. 5 non-woven-fabric layers and without exhalation valve. III category personal protective device (EN 149:2001/A1:2009).

HZ-KN95 - WHITE 5 PIECES

HZ-KN95 -63101374 - WHITE 1 PC

HZ-KN95-63101403 - BLACK 1PC

VISORS AND PROTECTION GLASSES



SURGICAL MASK WITH LACES AND VISOR - FPITV2R

Three-layer surgical face mask with laces and anti-fog visor. Made of TNT (non-woven fabric), it is seamless. CE marking in compliance with the European Directive 93/42 EEC (D.L. 46/97). Reference standard and certification EN ISO 14683: 2019 for type IIR masks.



RECLINING FACIAL PROTECTION SCREEN - L35101014

Personal protective equipment compliant with norm EN 166:2001. Reclining screen for face protection made of transparent polycarbonate. Adaptable to all head sizes thanks to a rear mechanical adjustment system. Effectively protects against liquids, particulates, and small mechanical impacts. Optic Class type 1. Thanks to the material with which it is made, the screen allows excellent visibility. The screen is reusable, and the front foam can be replaced by purchasing spare parts separately.



PROTECTIVE GLASSES - MED98704

Eye protection glasses. Transparent lens and frames. Scratchproof lens. Lateral protection on the bars. May be used



PROTECTIVE GLASSES - 511030100

Adjustable length and inclination of temples. Built-in side protectors. Frame made of hypoallergenic material. Protection against high-speed particles 45 m/s - low energy (F). Reference standards: EN166 EN170.

EYE WASHING KITS



EMERGENCY EYE WASHING KIT - CPS258

Recommended to all companies and on transport motor vehicles for emergency eye washing and its subsequent medication. Green case with wall support featuring: 1 empty eye washing bottle 500 ml; 1 eye washing sterile saline solution; wash 500 ml; 2 sterile eye pads; 2 sterile gauze packages 18x40cm; 1 pair of lister scissors 14,5cm; 1 plaster roll 1mx2cm; 1 towelette package; 1 hand mirror; 1 instruction manual.

Dimensions: about 25,5x9x19h cm.

Weight: about 1,2 kg.



EYE WASHING KIT

Eye washing kit featuring 2 disposable and sterile saline solution 500 ml bottles. Rapid eye pad. Thermoformed support with transparent frontal and powder free closing, equipped with small mirror, screws and nogs for wall mounting.

Dimensions: about 32x8x41 cm.

Weight: about 1,8 kg.

CPS202 – EYE WASHING KIT – EYE-SHAPED
TAP003 - EYE WASH CAP WITH ANTI-DUST COVER
SOL002 - PHYSIOLOGICAL SOLUTION 250ML - NOT INJECTION
SOL004 - PHYSIOLOGICAL SOLUTION 500ML - NOT INJECTION

DISINFECTANT PRODUCTS



UMONIUM MEDICAL WET WIPES 10 PCS - PF10748

Ultra-fast acting full-spectrum disinfectant, UMONIUM38 Medical Tissues are wet wipes which clean and disinfect, killing germs after 5 seconds of coming into contact with them. Cleaning and disinfection of non-invasive medical devices. The impregnating solution is validated for bactericidal, fungicidal, mycobactericidal (including B.K.), and viricidal action.



PHARMA STERIL DISINFECTING WIPES IN CANISTER - OX02

Disinfectant wipes in a canister, 100 pcs. Suitable for disinfection of surfaces and instruments. Contains quaternary ammonium salts. Effective against viruses: HIV, HBV, HCV and H1N1.



CHLOREXIDINE WIPE 10 PCS - TOW014

Disinfectant wipe based on chlorhexidine digluconate. Each sachet is produced with coupled film, specially made to contain the disinfectant liquid. Sachet dimensions: 60x80 mm. Wipe dimensions: 143x150 mm. Packaging: 10pcs.

DISINFECTANT PRODUCTS



PHARMA G DISINFECTANT

Disinfectant and detergent based on an association between didecyldimethylammonium chloride and orthophenylphenol. Can be nebulized with ULTRA LITE or in a one-liter bottle with trigger gun (available separately). Pleasant Lemon flavour. Surgical Medical Device.

Q641 - PHARMA G DISINFECTANT 1 LT
Q642 - PHARMA G DISINFECTANT 5 LT



SEPTALDEIDE FOR DISINFECTION AND STERILIZATION

Sterilizing solution containing glutaric surfactant aldehyde, with lemon essence. Efficient against spores and viruses including HIV, LAV, HBV and HCV. Pure application for sterilization of instruments and other objects, including thermolabile ones.

S431 - DISINFECTANT 1 LT
S432 - DISINFECTANT 5 LT



DISINFECTANT IMPREDIS SOFT - CE 19

Cleaning disinfectant for surfaces, surgical instruments, and dental impressions. Based on an association between two quaternary ammonium salts and chlorhexidine. Colorless and slightly orange scented. It does not contain alcohol and is therefore suitable for delicate surfaces and materials not compatible with alcohol. Efficient against HIV and HBV.

DISINFECTANT PRODUCTS



PHARMAIODIO DISINFECTANT

Povidone iodine 1000 ml bottle 10% P.M.C. iodine. Pharmaiodio is bactericide and virucidal and is suitable for pre-injection disinfection, pre-surgical preparation and intact skin antiseptic.

PF08110 - PHARMAIODIO 1 l
PF08111 - PHARMAIODIO 500 ml
PF08173 - PHARMAIODIO 125 ml



HYDROGEN PEROXIDE

1 lt. hydrogen peroxide bottle. Hydrogen peroxide, 3% solution – 10% volume. Stabilized with proper treatment. Suitable for skin cleaning.

MED98508 - HYDROGEN PEROXIDE 250 ml
MED98510 - HYDROGEN PEROXIDE 1 l



NON-INJECTABLE PHYSIOLOGICAL SOLUTION

0.9% non-injectable sodium chloride physiological solution. Sterile and apyrogenic solution in a recyclable polyethylene container. For wounds, suitable in-between and post-surgery. For bladder catheters and oral cavity cleaning. Cleaning and rinsing of surgical instruments and equipment.

195B230541 - PHYSIOLOGICAL SOL. 250 ml
195B230551 - PHYSIOLOGICAL SOLUTION 1 l



PHARMASIL 250 ML - B300

Rapid and persistent disinfectant for non-wounded skin. Mixture of didecyldimethylammonium chloride and orthophenyl. It rapidly gets rid of unwanted pathogenic microorganisms. Efficient against gram+ and gram- acid resistant bacteria. Scented, no burn, no stains. Surgical medical device.

DISINFECTANT PRODUCTS



SEPTAMAN GEL

Surgical Medical Device. Gel colourless alcohol solution (Ethyl alcohol 70%), ready to use, suitable for rapid and effective disinfection of intact skin. Kills bacteria, fungus, micro-bacteria and viruses. Furthermore, inactivates HIV, HBV, HCV and Rotaviruses. Contains no allergenic substances.

PF32459 - SEPTAMAN GEL 100 ml

PF32411 - SEPTAMAN GEL 500 ml

AN0425058 - WALL HOLDER SEPT. 500 ml



MEDICAL SPRAY DISINFECTANT - PF011

150 ml bactericide disinfectant. Bactericide and fungicide action, suitable for environmental disinfection treatments with automatic supply. This system allows the distribution of the product in the environmental thus efficiently disinfecting floors, walls, surfaces and air. It allows the complete and rapid disinfection of waiting, emergency, consulting and surgical rooms besides ambulances and classroom. Products diffusion is about 4 minutes long.



DISINFECTANT TOWELETTE - TOW223

Disinfectant towelette: each sachet contains 3 ml of Pharmaderm.

Dimensions of sachet: 65x95 mm circa.

Dimensions of towelette: about 140x190 mm.



4% AMMONIA WIPE - TOW013

Ammonia wipe in water solution 4%.
10 pcs package.

COLD/HOT THERAPY



DISPOSABLE INSTANT ICE PACK - FCT2000

Disposable ice with controlled temperature (0°C for about 30 minutes). Its immediate use during traumas avoids haematoma and swelling formation. It stops little haemorrhages and considerably reduces the pain. It can be kept without limits at ambient temperature.

Made in PVC material. Dimensions 17x13,5 cm. CE medical device.



DISPOSABLE NON-WOVEN ICE PACK - 02SF1853EASYICE

Disposable synthetic ice that keeps 0°C temperature for about 30 minutes.

Does not contain ammonium nitrate.

Made in non-woven material.

Dimensions: 13,5x17 cm.



SPRAY ICE

Spray ice bottle ready to use. Particularly suitable for sport activities, small traumas and contusions. CE medical device.

SP200DISPO24 - SPRAY 200 ml

SP400DISPO24 - SPRAY 400 ml

COLD/HOT THERAPY



HOT/COLD GEL - DISPOGEL

Reusable gelatin hot/cold Dispo GEL
It is a hot/cold reusable accumulator from outside source. Completely non-toxic, suitable for medical and domestic use: easily modeled, does not harden becoming ice and it is reusable. Cold therapy: keep in freezer for at least 2 hours. Take out of the freezer and put in a cloth and apply on the affected part. Hot therapy: place the bag in the microwave for at most 1 minute setting it to 650 Watt or dip the bag in boiling water for at most 5 minutes, remove and dry; in both cases insert the bag in a cloth and apply on the affected part.

Dimensions: about 14x24 cm.



DISPOSABLE INSTANT HOT PACK - 002SC1853HOT

Disposable instantaneous hot pack with a controlled temperature of 60°degrees lasting about 25 minutes.

Dimensions: about 13,5x17 cm.

ISOTHERMIC AND FIREPROOF BLANKETS

ISOTHERMIC BLANKET



Protection and rescue isothermic blanket
Prevents patients' hypothermia whether traumatized, burned or post-operated
Double function: warming and cooling
Silver side, if positioned toward the body, isolates and reflects, preserving about 80% of body warmth. Thus protects patient from hypothermia, coldness, freezing and shock.
Golden side: placed toward patient body and silver one, reflects sun ray maintaining body freshness. Non sterile.
Adult dimensions: about 160x210 cm.
Pediatric dimensions: about 100x80 cm.

COP216 - ADULT ISOTHERMIC BLANKET
FCT3010 - PEDIATRIC ISOTHERMIC BLANKET

BLIZZARD EMS BLANKET ORANGE - FCT3060



Designed specifically for the Emergency Services market the EMS is made from 2 layer Reflexcell™ and is widely used by all emergency services worldwide. The Blizzard EMS blanket works as an effective tool in the prevention and treatment of hypothermia. The Blizzard EMS blanket is a high performing blanket with self-adhesive closure that re-opens for easy access to casualties. The Blizzard EMS blanket is available vacuum packed for ease of storage: Small packaging. Easy to open and roll out. Size (vacuum packed): 12x11x5 cm
Weight: 250 g. Size (unpacked): 2,3x1,96 meters (fully stretched) 156.8cm (working)



FIREPROOF BLANKET

Fiberglass, fireproof blanket, EN 1869 compliant. Suitable for extinguishing small fires, especially for small oil fires.
Supplied with a case equipped with a ring for wall mounting and velcro opening.

COP208 - FIREPROOF BLANKET 120x180 cm
COP215 - FIREPROOF BLANKET 150x100 cm

ABSORBING TRANSFER SHEET WRASFER - 610007

TRANSFER SHEET

150 KG AND 3,5 LITERS

Wrasfer sheet in addition to being highly resistant (150 kg maximum carrying weight) allows the absorbency of fluids and exudates up to 3,5 liters, which guarantees the optimum containment of possible infections. Disposable device of the risk class I, made of 3-layer construction: PE-foil, nonwoven and cellulose fluff.

It has been specifically studied for the transfer assistance of patients and it has been designed for use in paramedic services, hospitals, hospital emergency wards, nursing homes, funeral chambers and for remains transportation in special conditions.



EASY TO USE

After its use the sheet can be easily eliminated and the carriers or stretchers are ready for the next use.

Sheet dimensions: 220x100 cm. Weight: 350g.

The application of the transfer sheet should be adequate to provide assistance to two/four patients.



WASTE SACKS AND SHEETS



VOMIT ABSORBING BAG - V60020

Vomit bag "100% Protected Hands."
Fitted with an absorbent pad with Alveol-Technology, it absorbs the vomit and allows its elimination: no dispersion. New MDPE liner: more resistant. Does not stick to the skin. Hypoallergenic.
Dimensions: 25x25 cm.
Packs containing: 20 pcs.
Latex Free, single use only.
CE Medical device.



ROLLED SHEETS - 005 P LU

Sheet rolls in pure cellulose, suitable for examination tables.
Embossed surface, two layers.
Average veil weight: 17 gr/mq.
Pack of 6 pieces.
Colour: white.
Dimensions: 70m x 60 cm.



WHITE VOMIT AND WASTE SACK - SAC052

White waste sack with hermetic closing. It takes little space and can be stored inside emergency bags and rescue rucksacks.
Dimensions: 30x50cm.
100 pieces package.



DISPOSABLE PILLOWCASE - FCT2466

Disposable non woven plasticized pillowcase, non-allergenic, reinforced, absorbing and impermeable.
Dimensions: 40x40 cm.
25 pieces package.

INSTRUMENTS - SCISSORS



FIRST AID SCISSORS CUT CLOTHING

Stainless steel first aid scissors cut clothing. Autoclavable plastic handle. Suitable for first aid bags and cases. Dimensions: 18 cm.

STR110010 - BLACK SCISSOR
STR110010/A - ORANGE SCISSOR
STR110010/B - BLUE SCISSOR
STR110010/G - YELLOW SCISSOR
STR110010/R - RED SCISSOR
STR110010/V - GREEN SCISSOR



LISTER SCISSORS

First aid cut-clothing stainless steel scissors. Suitable for emergency bags and cases .

STR110291 - LISTER SCISSORS 14 cm
STR110300 - LISTER SCISSORS 16 cm
STR110310 - LISTER SCISSORS 18 cm
STR110320 - LISTER SCISSORS 20 cm



LITTAUER LIGATURE SCISSOR - STR110330

Ligature scissors. Made in stainless steel
Dimensions: 14 cm.

INSTRUMENTS - SCISSORS



FIRST AID SCISSORS - FOR078

First aid scissors with plastic handle.
Stainless steel blades.
Dimensions: 10 cm.



STRAIGHT SCISSORS WITH BLUNT

Straight medication scissors with blunt tips.
Made in stainless steel.

STR110080 - SIZE 14 cm
STR110090 - SIZE 16 cm



STRAIGHT SCISSORS WITH SHARP TIPS

Straight medication scissors with sharp tips.
Made in stainless steel.

STR110040 - SIZE 14 cm
STR110050 - SIZE 16 cm



STRAIGHT SCISSORS WITH ALTERNATE TIPS

Straight medication scissors with alternate tips, blunt/sharp. Made in stainless steel.


STR110060 - SIZE 14 cm
STR110070 - SIZE 16 cm



SCISSORS WITH SHARP TIPS

Medication scissors with bent alternate tips, blunt/sharp. Made in stainless steel.

STR110030 - SIZE 14 cm
STR110032 - SIZE 16 cm

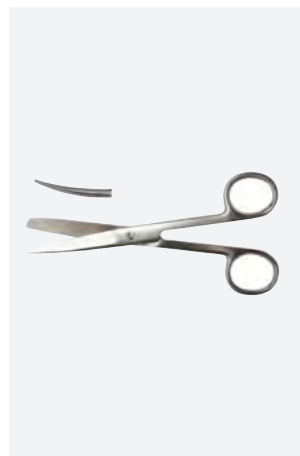

BENT SCISSORS WITH ALTERNATE TIPS

Medication scissors with bent alternate tips, blunt/sharp. Made in stainless steel.

STR110020 - SIZE 14 cm
STR110035 - SIZE 16 cm

INSTRUMENTS - SCISSORS



MULTIPURPOSE EMERGENCY SCISSORS - FOR299

Multipurpose Emergency scissors
Ideal in cases of emergency to break crystals, cut clothes, cut safety belts, cut thick objects (knife function). Made of metal that is not subject to oxidation, without plastic parts. You can open it completely, dividing it into two pieces, to wash it more easily. Completely autoclavable. Supplied with cover and accessories for wall mounting.



ROBIN SCISSORS - 124900

Multi-function emergency scissors, Robin original. Made in Solingen steel. It easily cuts fasten belts. Its handle allows car glasses crush. Supplied with wall support. Weight: 400 g.

124900 - ROBIN SCISSORS

124901 - WALL HOLDER FOR ROBIN SCISSOR



STERILE DISPOSABLE SCALPEL

Disposable sterile scalpel, singularly packaged. Handle in abs material. Stainless steel blade. Blade protection in polythene, easy to remove.

90466 - DISPOSABLE SCALPEL SIZE 10

90467 - DISPOSABLE SCALPEL SIZE 11

90468 - DISPOSABLE SCALPEL SIZE 12

INSTRUMENTS



REMOVAL PENCIL - MAT022

Small metallic splinters eye tool remover. One side you have a NEODIMIO magnetic point, which enable you to withdraw all metallic splinters avoiding any eventual dangerous eye approach. On the other side you have an non-abrasive nylon loop which is used for non metallic splinters. Main body and support made of anodised aluminium and the neoprene seals guaranty a better pencil hold. Supplied in a practical box. CE medical device.



PERCUSSION HAMMER - STR110270

Percussion hammer. Zinc plated handle. Head in pvc. Complete with needle. Dimensions: about 18x6x1,5 cm.



DISPOSABLE TRICHOTOMY RAZOR - RSA225

Trichotomy razor suitable for hairs removal in the area of the skin supposed to surgical operation so to avoid possible infections. Stainless steel blade. Disposable non-sterile single blade.

INSTRUMENTS - FORCEPS



MAGILL FORCEPS


Intubation Magill forceps. Made in autoclavable and disinfectable stainless steel.

STR110230 - MAGILL FORCEPS 16 cm
STR110240 - MAGILL FORCEPS 20 cm
STR110250 - MAGILL FORCEPS 25 cm



CLAMP PLASTIC FORCEPS - STR110450

Clamp plastic forceps.
Dimensions: 12 cm.

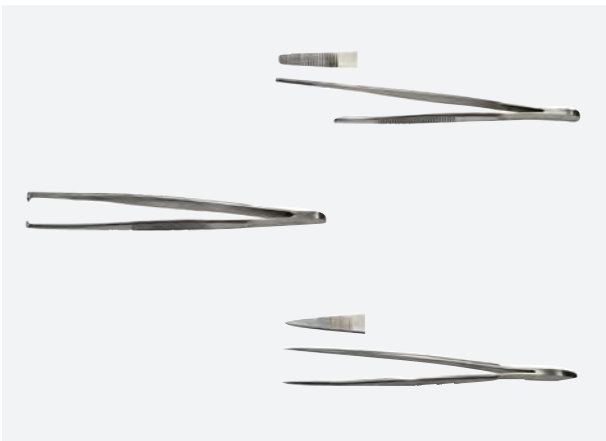

SURGICAL AND ANATOMIC TWEEZERS

Anatomic tweezer made in stainless steel.

STR110100 - ANATOMIC TWEEZER THIN 14 cm
STR110110 - ANATOMIC TWEEZER 14 cm
STR110120 - ANATOMIC TWEEZER 16 cm
STR110130 - SURGICAL TWEEZER 14 cm
STR110140 - SURGICAL TWEEZER 16 cm



PEAN, KOCHER AND WIRE CONDUCTOR FORCEPS

Tweezer made in stainless steel.

STR110150 - PEAN FORCEPS 14 cm
STR110160 - PEAN FORCEPS 16 cm
STR110170 - KOCHER FORCEPS 14 cm
STR110180 - KOCHER FORCEPS 16 cm
STR110190 - MATHIEU WIRE CONDUCTOR
FORCEPS 14 cm
STR110200 - MATHIEU WIRE CONDUCTOR
FORCEPS 16 cm
STR110220 - WIRE CONDUCTOR FORCEPS
14 cm
STR110280 - MOSQUITO FORCEPS

INSTRUMENTS - FORCEPS



RING CUTTER FORCEPS - STR110500

Stainless steel ring cutter forceps.
Dimensions: 18 cm.



NYLON FORCEPS - PIN 103

Nylon forcep singularly packed
Sterile. 8 cm.



STAINLESS STEEL TONGUE FORCEPS - STR110210

Stainless steel tongue forceps.
Dimensions: 17 cm.



PLASTIC TONGUE FORCEPS - PIN201

Polycarbonate tongue forceps.
Dimensions: 16 cm.



PLASTIC TICK-REMOVING TWEEZERS - STR110462

Suitable for animals. The pen-shaped structure allows to capture the body of the tick in order to extract it with a rotatory movement. The correct extraction prevents infections.
Dimensions: 9 cm.
Colours availability: random.

PROTECTION KITS



BASIC PROTECTION KIT - VAL44094

Protection kit including:
20 x FFP2 protection masks
1 x 100 ml disinfectant hand gel.



PROTECTION KIT PLUS - VAL44091

Protection kit in pouch including:
10 x FFP2 protection masks
1 x finger pulse oximeter
1 x digital infrared thermometer
1 x 500 ml disinfectant hand gel
1 pack of 100 disinfectant tissues for surfaces and medical devices
Dimensions: approx. 25x15x12 cm



PROTECTION KIT WITH THERMOMETER - VAL44092

Protection kit in pouch including:
20 x FFP2 protection masks
1 x digital infrared thermometer
1 x 500 ml disinfectant hand gel
1 pack of 100 disinfectant tissues for surfaces and medical devices
Dimensions: approx. 25x15x12 cm



PROTECTION KIT WITH PULSE OXIMETER - VAL44093

Protection kit in pouch including:
20 x FFP2 protection masks
1 x finger pulse oximeter
1 x 500 ml disinfectant hand gel
1 pack of 100 disinfectant tissues for surfaces and medical devices
Dimensions: approx. 25x15x12 cm

EMERGENCY KITS



PERIPHERAL VEINS CANNULATION SET - VAL44060

Peripheral veins cannulation set in pouches with transparent frontal, containing:
1 haemostatic belt, 5 canula needles (16/17/18/20/22), 3 butterfly needles (21/23/25), 1 silky plaster roll 2,5 cm, 4 disinfectant wipes, 2 IV administration sets.
Dimensions : about 27x12x9(h) cm



VALVE KIT TYPE HEIMLICH - 220000066910

Non-return valve for thoracic drainage (Heimlich type). Valve body in flexible and transparent polyurethane with natural rubber sleeve in “duck’s bill” shape acting as a non-return valve. This device is connected directly onto the patient’s pleural or thoracic drain (adult, paediatric), or indirectly, using connector code 800.01 and an extension tube for neonates and infants. This product allows thoracic drainage to continue without risk to the patient even when transported, as the non-return valve means that the drain does not have to be clamped.



PLEUROTHORAX SET - 250.200

Closed System comprising: secretion bag 2,0 ltr. with non-return valve, stopcock on the base and three-way stopcock (fixed). Transfusion Tube Luer-Lock (unscrewable). Puncture Cannula Diameter 2 mm; puncture length 75 mm. 50-ml-syringe Luer-Lock. Sterile double wrapped; Peel-Pack.

FIRST AID KITS



TRAVEL KIT - VAL44081

Travel first aid kit complete with: 2 small bags for minor burns treatment, 2 2-pieces non-woven sterile compress, 1 plaster strip, 2 glove pairs in non-sterile package, 3 disinfectant tissues, 1 nylon non sterile tweezers, 1 cohesive elastic bandage, 1 first aid scissors, 2 ammonia towels, 1 suture plasters package, 1 disposable standard instant ice, 40 plasters.

Dimensions: about 15x12x8,5 cm.



PET KIT - VAL44076

First aid kit ideal during small emergencies involving pets, designed especially for dogs. Small handbag in washable PVC made of coloured waterproof and tearproof fabric with side Velcro strip for anchoring to backpacks or rescue bags. The handbag contains the following material: Plaster roll. Cohesive elastic dressing bandage. DIN medication package. 2x sterile gauze (2 pcs). 2x pairs of gloves. 2x disinfectant wipes. Small first aid scissors. Isothermal emergency blanket. Sterile single-use syringe. Nylon tweezers. Tick-removing tweezers. Dimensions: about 10x7x13(h) cm.



MOUNTAIN RESCUE FIRST AID KIT - VAL44082

Mountain rescue first aid kit complete with: 1 emergency isothermal blanket, 2 small bags for minor burns treatment, 5 2-pieces non woven compresses, 1 sterile sheet, 1 bag of aluminized dressing, 1 striped medicated plaster, 2 pairs of non-sterile gloves, 4 disinfecting wipes, 1 small first aid scissors, 1 sterile nylon tweezers, 1 life saver wipe, 2 napkins, ammonia 4%, 1 pack of suture plasters, 1 pack of standard disposable ice, 40 assorted plasters, 1 instant hot compress, 250ml of non-injectable saline solution. Dimensions: about 23x14x11cm.

FIRST AID KITS



DIN13164 PACKAGE - CPS685

Red nylon handbag. Content complies with what stated by norm DIN 13164 (mandatory on automotive vehicles in certain European Nations). 8 adhesive bandages 10x6cm. Plaster roll 5 m x 2,5 cm. 2 gauze bandages 6 cm x 4 m. 3 gauze bandages 8 cm x 4 m. 1 pair of bandage-cutting scissors. Package of 4 vinyl gloves. Non-woven fabric sheet 60x80 cm. 2 non-woven fabric sterile sheets 60x40 cm. 3 sterile swab bandages 80x100 mm. Sterile swab bandage 100x120 mm. 6 non-woven fabric sterile swabs 10x10 cm. 2 non-woven fabric triangular sheets 96x96x136 cm. Isothermal blanket. Dimensions: 25x15x7,5 cm. Weight: 0,45 kg.



DEFI SET - MED98950

Ready to use set for early defibrillation, consisting of:

- Clothes-cutting scissors
- Trichotomy razor
- 1 pair of disposable gloves
- 2 disinfectant wipes
- 2 non-woven absorbing compresses 10x10 cm
- Life key ventilation mask

FIRST AID KITS



FIRST AID INDIVIDUAL KIT - VAL44078

Complete first aid and medication kit:
Isothermic blanket, 3,5 ml Idrogel sachet,
2 sterile non woven compresses, Sterile
medication cloth 60x80 cm, 2 Alumed
medication pack size M, Transparent plaster
roll 2,5 cm, Cotton medicated plaster 10x6
cm, 2 pairs of gloves in non-sterile package,
2 disinfectant wipes, Small first aid scissors,
Nylon sterile tweezer, Life saver, Ammonia
4% wipe, First aid instructions.
Dimensions: about 15x12x8,5 cm.
Weight: about 200 g.



CAR DIN 13164 SET INCREASED - VAL44362

Practical first aid bag for cars: 1 plaster
roll, 8 strip plaster, 1 wound medication
package L, 3 wound medication packages
M, 2 immobilization triangular clothes, 1
sterile medication cloth for burn, 2 sterile
medication clothes for burns, 3 big sized
medication bandage gauzes, 2 medium
sized bandages, 6 sterile gauze compresses,
1 angular medication scissors, 1 emergency
isothermic blanket, 4 vinyl protection gloves,
8 disinfectant wipes, 3,5 ml Idrogel burns
sachet, 2 4% ammonia wipes for insect
stings.
Dimensions: about 27x12x9(h) cm.
Weight : about 0,5 kg.



FIRST AID KIT-MOTORBIKE DIN 13167 - 17008

First Aid Kit for motorbikes with content
according to DIN 13167. 1 roll of adhesive
plaster, 14-piece adhesive dressing
assortment: 4 wound dressings, 2 fingertip
dressings, 2 finger bandages, 2 plasters, 4
plasters. 1 first aid dressing, sterile, 1 first
aid dressing, sterile, 1 burn dressing, sterile,
1 rescue blanket, 1 bandage scissors, 2
pairs of disposable gloves, 2 wet cleansing
wipes for unwounded skin, 1 first aid
instructions, 1 list of contents.

REPLANTATION KITS



REPLANTATION SET - VAL40954

Replantation set in a tearproof impermeable equipped bag with 4 coloured pouches with handles transparent frontal, complete of:

1 sterile replacement/transport set for finger,
1 sterile replacement/transport set for hand,
1 sterile replacement/transport sack for arm,
1 sterile replacement/transport sack for leg,
Hemostatic belt/lace, 1 tourniquet type Martin,
10 sterile gauze packs 10x10cm,
1 plaster roll 2,5cm, 1 emergency pad bandage,
1 Alumed cloth 40x60cm, 1 anatomic forceps,
1 surgery forceps, 1 Pean forceps.

Dimensions: about 48x16x26(h) cm.

Weight: about 5,5 kg.

MED98820 - HAND REPLANTATION SET - 300x250 mm
MED98824 - ARM REPLANTATION SET - 700x250 mm
MED98828 - LEG REPLANTATION SET - 1100x400 mm
MED98832 - FINGER REPLANTATION SET - 200x200 mm
MED98810 - HAND REPLANTATION BAG - 300x200 mm
MED98811 - FINGER REPLANTATION BAG - 200x200 mm
MED98812 - ARM REPLANTATION BAG - 700x200 mm
MED98813 - LEG REPLANTATION BAG 1100x400 mm
VAL40954 - REPLANTATION SET

REPLANTATION LIMBS SET CASE - VAL40950

Replantation limbs set case, made in shockproof orange ABS, featuring:

1 sterile replacement/transport set for finger,
1 sterile replacement/transport set for hand,
1 sterile replacement/transport sack for arm,
1 sterile replacement/transport sack for leg,
1 haemostatic lace/belt, 10 sterile gauze packs 10x10cm,
1 plaster roll 2,5cm, 1 emergency pad bandage,
1 Alumed cloth 40x60cm, 1 anatomic forceps,
1 surgery forceps, 1 Pean forceps.

Dimensions: about 42x15x31(h) cm.

Weight: about 7,5 kg.



EMERGENCY DELIVERY KITS



EMERGENCY DELIVERY SET - VAL40730

Emergency delivery set in pouch with transparent frontal, featuring: 1 clamp forceps, 2 sterile umbilical forceps, 3 sterile medication bandages, 1 mucus aspirator, 1 baby isothermic blanket, 2 sterile medication clothes 40x60cm, 2 sterile highly absorbing compresses, 5 sterile gauze compresses 10x10cm, 2 elastic bandages, 1 plaster roll, 1 pair of sterile gloves, 1 first aid scissors, 1 bottle of Povidone Iodine.
Dimensions: about 27x12x9h cm.



EMERGENCY DELIVERY SET IN CASE - VAL40732

Emergency birth set in case with transparent frontal, featuring: 1 Clamp forceps, 2 sterile umbilical forceps, 3 sterile medication bandages, 1 mucus aspirator, 1 baby isothermal blanket, 2 sterile medication clothes 40x60cm, 2 sterile highly absorbing compresses, 5 sterile gauze compresses 10x10cm, 2 elastic bandages, 1 plaster roll 1 pair of sterile gloves, 1 scissors 14cm, 1 bottle of Povidone Iodine.

Dimensions: about 25,5x9x19,5(h) cm.

WATER-JEL BURN DRESSING

Water-Jel Burn Dressings provide fast, one-step relief for emergency first aid when someone has been burned. The dressings stop burn progression and cool the skin, relieve pain and protect against airborne contamination. Water-Jel Burn Dressings are made with a scientifically formulated cooling gel and utilize a medical-grade nonwoven pad that is non-adherent, so they will not stick to burn injuries. The dressings are soaked in a water-based, water-soluble gel and sterilized to ensure a contamination-free atmosphere for emergency treatment of burns. With Water-Jel Burn Dressings, it takes no special training to give fast and effective emergency burn care (4-6 hours) for any type of burn injury. All Water-Jel Burn Dressings have a 5-year shelf life from date of manufacture.



WATER-JEL BURN DRESSING

B0206 - BURN DRESSING 5x15 cm
B0404 - BURN DRESSING 10x10 cm
B0416 - BURN DRESSING 10x40 cm
B0818 - BURN DRESSING 20x45 cm
B0820 - HAND BURN DRESSING 20x55 cm
B1216 - FACE BURN DRESSING 30x40 cm



WATER-JEL IN GEL

GBJ140 - BURN JEL 4.0 gr
GBJ80 - BURN JEL 80 ml



WATER-JEL FIRE BLANKET

Fire blanket Water-Jel provide fast, one-step relief for emergency first aid when someone has been burned. The dressings stop burn progression and cool the skin, relieve pain and protect against airborne contamination. It takes no special training to give fast and effective emergency burn care. It can be used to extinguish flames and smoldering clothes.

G3630P - FIRE BLANKET 91x76 cm
G7260C - TUBE FIRE BLANKET 183x152 cm

BURN DRESSING KITS



MINI BURN KIT - VAL 40516

Mini burn kit for burns treatment contained in a pouch, featuring: 2 WATER-JEL burn medication gauze 10x10 cm. 1 WATER-JEL burn medication gauze 20x20 cm. 2 small towelettes for minor burns, 3.5 ml each. 1 pair of sterile gloves. 1 elastic medication bandage. 1 emergency scissors with plastic handles. Dimensions: about 23x13x10 cm ca. Weight: 0.7 kg ca.



IDROGEL BURN PROOF SPRAY NO GAS 125 ML - 500520

125 ml gel spray no gas bottle. Brustop burn proof provides the necessary skin hydration, it cools and helps burn depth decrease, rash reduction and immediate pain release. IIb risk class device, clinically tested, suitable for use on I,II and III level scalding even on wounded skin. Different sizes and dimensions, moldable to all needs. Sterile packaging. Non-toxic and causes no irritation. Suitable for pediatric use. Does not stick to the skin, easily removable, elastic and flexible. Biodegradable.


ITALY MEDICATION CASE D.M. 388/2003 - D.LGS. 81/2008 ANNEX 1



ITALY MEDICATION CASE ANNEX 1

Medication case complete with:
5 pairs of disposable sterile gloves; 1 anti-squirt visor; 1 litre of cutaneous Povidone-iodine with 10% iodine; 1.5 litre of saline solution; 10 packages of sterile gauzes 10 x 10 cm – single packaging; 2 packages of sterile gauzes 18 x 40 cm – single packaging; 2 disposable sterile sheets; 2 disposable sterile medical forceps; 1 medium-sized elastic bandage; 1 package of hydrophilic cotton; 2 packages of various sized plasters ready-to-use; 2 rolls of high plaster 2,5cm; 1 pair scissors; 3 tourniquets; 2 packages of instant ice; 2 disposable sanitary waste bags; 1 thermometer; 1 sphygmomanometer; 1 stethoscope.
Dimensions: about 39,5x13,5x26(h) cm.
Weight: about 5,8 kg.

VAL45180 - MEDICATION CASE D.L. 81/2008 ANN. 1
VAL45184 - MEDICATION CABINET D.L. 81/2008 ANN. 1
VAL45150 - MEDICATION BAG 33 D.L. 81/2008 ANN. 1
MED98878 - MEDICATION PACKAGE D.L.81/2008 ANN. 1 BASE
MED98874 - MEDICATION PACKAGE D.L.81/2008 ANN. 1 COMPLETE

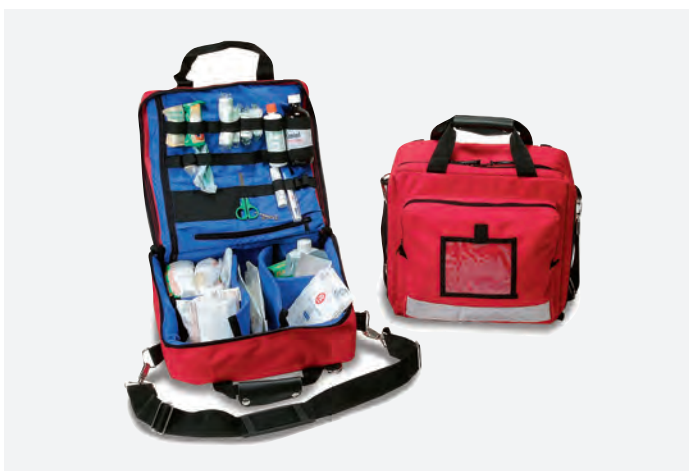

D.LGS. 81/2008 ANNEX 2



ITALY MEDICATION CASE ANNEX 2

Medication case complete with:
2 pairs of disposable sterile gloves; 125 ml of cutaneous Povidone Iodine with 10% iodine; 250 ml of saline solution; 3 packages of sterile gauzes 10 x 10 cm – single packaging; 1 packages of sterile gauzes 18 x 40 cm – single packaging; 1 disposable sterile medical forceps; 1 package of hydrophilic cotton; 1 package of various sized plasters ready-to-use; 1 roll of high plaster 2,5 cm; 1 roll of elastic edged bandage – height 10 cm; 1 pair scissors; 1 tourniquet; 1 package of instant ice; 1 disposable sanitary waste bags;
Dimensions: 25x9x19,5 C cm.
Weight : about 1,3 kg.

VAL45170 - MEDICATION CASE D.L. 81/2008 ANN. 2
VAL45174 - MADICATION POUCHE D.L. 81/2008 ANN. 2
MED98870 - MEDICATION PACKAGE D.L.81/2008 ANN.2



BOLZANO MEDICATION CASE D.P.P. 25/2005 ANNEX 1

BOLZANO MEDICATION CASE ANNEX 1



Medication case complete with:
5 pairs of disposable sterile gloves; 1 anti-squirt visor; 1 litre of cutaneous Povidone Iodine with 10% iodine; 1.5 litre of saline solution; 10 packages of sterile gauzes 10 x 10 cm – single packaging; 2 packages of sterile gauzes 18 x 40 cm – single packaging; 2 disposable sterile sheets; 2 disposable sterile medical forceps; 1 medium-sized elastic bandage; 1 package of hydrophilic cotton; 2 packages of various sized plasters ready-to-use; 2 rolls of high plaster 2,5cm; 1 pair scissors; 3 tourniquets; 2 packages of instant ice; 2 disposable sanitary waste bags; 1 thermometer; 1 sphygmomanometer; 3 medication packages - 8 cm; 3 medication packages - 10 cm; 3 elastic edged bandage - size 8 cm; 3 elastic edged bandage - size 10 cm; 1 isothermic blanket; 1 protection handkerchief for artificial ventilation; 250 ml bottle of disinfectant; 1 stethoscope.
Dimension: about 39,5x13,5x26 cm.
Weight: about 6,2 kg.

VAL45186 - BZ MEDICATION CASE D.L. 81/2008 ANN. 1
VAL45188 - BZ MEDICATION CABINET D.L.81/2008 ANN. 1
VAL45152 - BZ MEDICATION BAG 33 D.L. 81/2008 ANN. 1
MED98868 - BZ MEDICATION PACKAGE D.L.81/2008 BASE ANN. 1
MED98866 - BZ MEDICATION PACKAGE D.L.81/2008 COMPLETE ANN. 1



BOLZANO MEDICATION CASE D.P.P. 25/2005 ANNEX 2

BOLZANO MEDICATION CASE ANNEX 2



Medication case complete with:
2 pairs of disposable sterile gloves; 125 ml of cutaneous Povidone Iodine with 10% iodine; 250 ml of saline solution; 3 packages of sterile gauzes 10 x 10 cm – single packaging; 1 packages of sterile gauzes 18 x 40 cm – single packaging; 1 disposable sterile medical forceps; 1 package of hydrophilic cotton; 1 package of various sized plasters ready-to-use; 1 roll of high plaster 2,5 cm; 1 roll of elastic edged bandage – height 10 cm; 1 pair scissors; 1 tourniquet; 1 package of instant ice; 1 disposable sanitary waste bags; 2 medication packages - 8 cm; 2 medication packages - 10 cm; 2 elastic edged bandage - size 8 cm; 2 elastic edged bandage - size 10 cm; 1 isothermic blanket; 1 protection handkerchief for artificial ventilation.
Dimensions: about 25,5x9x19,5 cm.
Weight: about 1,5 kg.

VAL45176 - BZ MEDICATION CASE D.L. 81/2008 ANN.2
VAL45178 - BZ MEDICATION POUCHE D.L. 81/2008 ANN.2
MED98864 - BZ MEDICATION PACKAGE D.L.81/2008 ANN.2





MEDICAL OXYGEN

OXYGEN THERAPY MASKS



OXYGEN THERAPY MASK- 11914020304

Adult oxygen therapy mask. Manufactured in transparent, atoxic and soft material complete with adjustable elastic. Supplied with 210 cm crushproof tube . Singulary packaged.



PEDIATRIC OXYGEN THERAPY MASK - 11914020302

Pediatric oxygen therapy mask made of transparent, atoxic and soft material complete with adjustable elastic. Supplied with 210 cm crushproof tube. Singulary packaged.



ADULT OXYGEN THERAPY MASK WITH RESERVOIR - 11914020313

Oxygen therapy mask with reservoir and non-return valve. Made of transparent, atoxic and soft material complete with adjustable elastic. High concentration (60% Oxygen). Metal nose clip allowing a safe and comfortable appliance. Supplied with 210 cm crushproof tube and universal connector. Size: adult nr. 6.



PEDIATRIC OXYGEN THERAPY MASK WITH RESERVOIR - 11914020311


Oxygen therapy mask with reservoir and non-return valve. Made of transparent, atoxic and soft material complete with adjustable elastic. High concentration (60% Oxygen). Metal nose clip allowing a safe and comfortable appliance. Supplied with 210 cm crushproof tube and universal connector. Singulary packaged.

OXYGEN THERAPY MASKS



ADULT VARIABLE CONCENTRATION MASK- 11914020328

Mask for oxygen therapy with variable concentration. Made of transparent, non-toxic and soft material, complete with adjustable elastic. Metal nose clip for safe and comfortable application. Equipped with a 210 cm anti-crushing hose and an adjustable concentration connection. Size: adults.

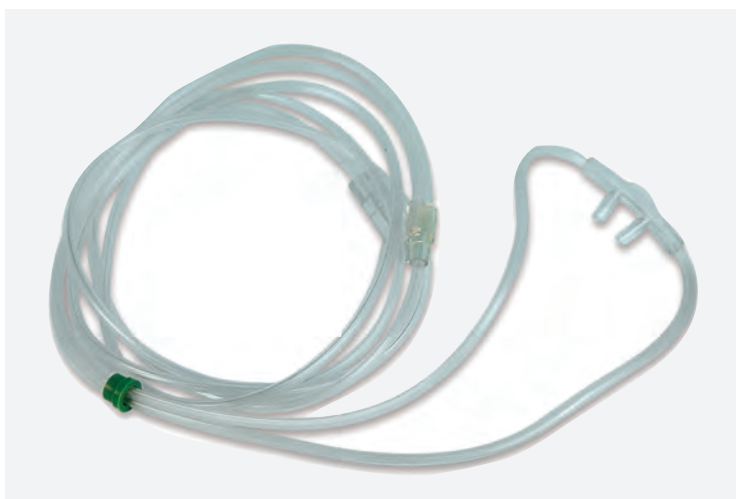

OXYGEN THERAPY NASAL CANULA 180CM - 11914020101

Disposable oxygen therapy nasal cannula. Manufactured in atoxic and soft material. It allows freedom of movement and do not interfere with feeding and vocal activities. Supplied with a 180 cm tube in singularly sealed package.



CRUSHPROOF OXYGEN FOR RESUSCITATOR BAGS - 11914020204


Oxygen 2 m crushproof tube extension complete with male/female connectors.

OXYGEN DISTRIBUTION KITS FOR AMBULANCES



Our medical oxygen distribution kit has been carefully designed and built for emergency vehicles, more specifically ambulances. The modularity of the system allows you to choose the number and type of terminal connections required, and the inclusion or not of the pressure gauge and/or of the manual oxygen cylinder switch. In fact, we offer kits with specific codes, depending on the different configurations (terminal connections, length of the tubes and tubes connections). We do also offer a plate that will allow you to fit the oxygen distribution kit under the inner wall of the ambulance.

Oscar Boscarol Srl and its official service centres are the only ones authorized for the repair, maintenance and technical inspection of Boscarol devices: contact us for information.

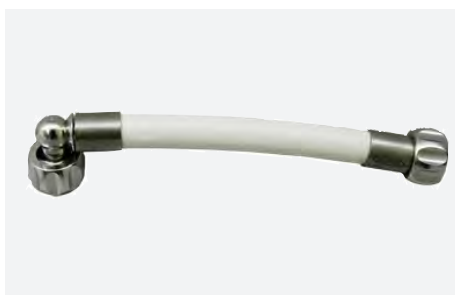

SWITCH



EXTERNAL SWITCH

Mechanical manual switcher for oxygen control on oxygen distribution kits with 2 inlets and 1 output. Both available with lateral and back inlets/outlets. It allows 2 low pressure oxygen flow selections (max. 8 bar – 800 kPa). Oxygen compatible, made with properly treated material on surface so to grant full compliance to ISO15001. Class of risk: IIa. Big control knob easy to use also with gloves. 3/8" oxygen hose connector. Optional: frontal plate for wall embedding.

OXIP7000 - EXT SWITCH 2EL-1UL
OXIP7001 - EXT SWITCH 2EP-1UP
OXIP7003 - EXT SWITCH 2EP-1UL
OXIP7004 - EXT SWITCH 2EL-1UP



EXTERNAL SWITCH WITH FRONTAL PLATE

Mechanical manual switcher for oxygen control on oxygen distribution kits with 2 inlets and 1 output. Both available with lateral and back inlets/outlets. It allows 2 low pressure oxygen flow selections (max. 8 bar – 800 kPa). Oxygen compatible, made with properly treated material on surface so to grant full compliance to ISO15001. Class of risk: IIa. Big control knob easy to use also with gloves. 3/8" oxygen hose connector. Frontal plate for wall embedding.

OXIP7000P - EXT SWITCH 2EL-1UL PLATE
OXIP7001P - EXT SWITCH 2EP-1UP PLATE
OXIP7003P - EXT SWITCH 2EP-1UL PLATE
OXIP7004P - EXT SWITCH 2EL-1UP PLATE

FLOWMETER AND HUMIDIFIER

AUTOCLAVABLE HUMIDIFIER TR200 9/16" 200ML - 000070100

Reusable humidifier for wall and plug-in flowmeters. Increases the percentage of humidity and promotes better patient breathing. Made of polycarbonate frame and chrome-plated brass fittings. Autoclavable (121°C 15 min.) Internal volume: 200 ml. Maximum applicable pressure value 500 kPa. Maximum applicable flow value 20Lpm. Dimensions: 78x60x192(h) mm. Weight: 0.138 kg. Humidification capacity with oxygen supply (relative humidity at 18.7°C): 14% inlet, 89% outlet, 73% after 2 m of hose with Ø5 with mask.

AUTOCLAVABLE HUMIDIFIER 9/16 - 300 ML - 000075010

Disposable humidifier for wall-mounted or plug-in flowmeters. Increases moisture content and promotes better patient breathing. It is a disposable device. Made with ABS frame and lightweight polycarbonate jar. Internal volume: 300 ml. Maximum applicable pressure value 500 kPa. Maximum applicable flow value 20Lpm. Humidification power with oxygen supply (relative humidity at 18.7°C): 14% inlet, 89% outlet, 73% after 2 m of hose with Ø5 with mask. Dimensions: 83x70x178(h) mm. Weight: 0.076 kg

FLOWMETER 0-15 LPM

OXI0566 - FLOWMETER 0-15 LPM UNI INPUT
OXI0567 - FLOWMETER 0-15 LPM AFNOR INPUT
OXI0568 - FLOWMETER 0-15 LPM DIN INPUT



PRESSURE REDUCER

PRESSURE REDUCER P300 UNI - 4501051

Pressure reducer for medical oxygen compliant with the standard ISO 10524-1
Made of surface treated brass
P1 (maximum inlet pressure): 300 bar
P2 (nominal output pressure): 4,5 bar
Reducer capacity at 200 bar (20.000 kPa): max. 250 lpm (litres per minute)
Oxygen cylinder connection type UNI 4406
Output connection: 3/8" (three eighths of an inch)
Compliance: CE marking according to MDD 93/42/CEE
Oxygen compatible according to ISO 15001
Operating temperature: -20÷60 °C
Storage and conservation temperature: -20÷60 °C
Warranty: 2 years from date of purchase (see warranty conditions)
Periodic safety check: 5 years from the date of manufacture
Dimensions: about 110x51x58mm
Weight: about 580g





MEDICATION SYMBOL - STE 141510

Adhesive reflex symbol.
Dimensions: 4x4 cm.



INSTRUMENTS SYMBOL - STE 141512

Adhesive reflex symbol.
Dimensions: 4x4 cm.



SUCTION UNIT SYMB. - STE 141514

Adhesive reflex symbol.
Dimensions: 4x4 cm.



OXYGEN BOTTLE SYMB. - STE 141516

Adhesive reflex symbol.
Dimensions: 4x4 cm.



LARYNGOSC. SYMBOL - STE 141518

Adhesive reflex symbol.
Dimensions: 4x4 cm.




RESUSCITAT. SYMBOL - STE 141520

Adhesive reflex symbol.
Dimensions: 4x4 cm.



PHIALS SYMBOL - STE 141522

Adhesive reflex symbol.
Dimensions: 4x4 cm.

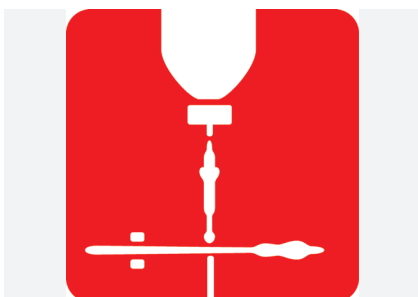

INFUSIONS SYMBOL - STE 141524

Adhesive reflex symbol.
Dimensions: 4x4 cm.



GLOVES SYMBOL - STE 141526

Adhesive reflex symbol.
Dimensions: 4x4 cm.



DIAGNOSTIC SYMBOL - STE 141528

Adhesive reflex symbol.
Dimensions: 4x4 cm.



INTERNATIONAL SYMBOLS - STE 141500

Adhesive reflex international
symbolisms. 12 different pieces
packages. Dimensions: 4x4 cm.

COVERALL AND SHOE COVER

HEAT-SEALED SEAM PROTECTIVE COVERALL 4B 5B 6B-XL - 21083



White disposable coveralls with blue heat-sealed seams, type 4B, 5B, 6B. Applications: ambulatory services of medical institutions, wards, laboratories, and other uses for regular insulation. Breathable microporous non-woven fabric coverall, with hood, heat-sealed seams, elasticized cuffs, elasticized waist, front zip. Dustproof, resistant to dry particles/light splashes/ viruses and blood permeate. TYPE-4B: splash-proof (spray), offers protection against liquid chemical saturation. TYPE-5B: dust-tight, offering full-body protection against airborne solid particles. TYPE-6B: limited liquid splash-proof, offering limited protection against light liquid splashes. Category III PPE, compliant with Regulation (EU) 2016/425. Certified according to EN 14126:2003, EN 1149-5:2008, EN 1073-2:2002 excluding 4.2, EN 14605:2005/A1:2009, EN ISO 13982-1:2004/A1:2010; EN 13034:2005/A1:2009. Dimensions: size XL, chest approx. 108-116 cm, height approx. 182-188 cm. Weight: 60 g approx. Individually packaged.

PERFECTO PE SHOE COVER 100 PCS - CCPE1200PER



Blue disposable shoe cover in PE (polyethylene), available in one-size, with ankle elastic. Thin, lightweight and waterproof. Ideal for both patients and operators. Medical Device Class I according to European Directive 93/42/EEC and following amendments. Storage: room temperature in a dry place, away from heat sources, flames, or sparks. Packaging: box of 100 pcs.

OTHERS



FACE PROTECTIVE SHEETS - DID150198

Ventilation's handkerchief with filter protection. It allows the mouth to mouth ventilation in absolute safety. Transparent PVC sheet, single-use, non-toxic. Suitable for training on manikins.

Dimensions: about 19,5x25 cm.
Package of 50 pieces in roll.




TRIAGE TAGS - PC190000

Triage cards made of cardboard. The triage card is composed not only from the standard part reporting the state and the practiced therapy but also by 4 stubs, perforated tear-off, with 4 different colors, indicating the severity of the patient. Package of 25 pieces
Dimensions: 10,5x20,5 cm.



SAFETY SIGN OF FIRST AID - ACC170550


Safety Sign of First Aid (white cross on a green background).
Made of aluminum.
Dimensions: 25x25 cm.

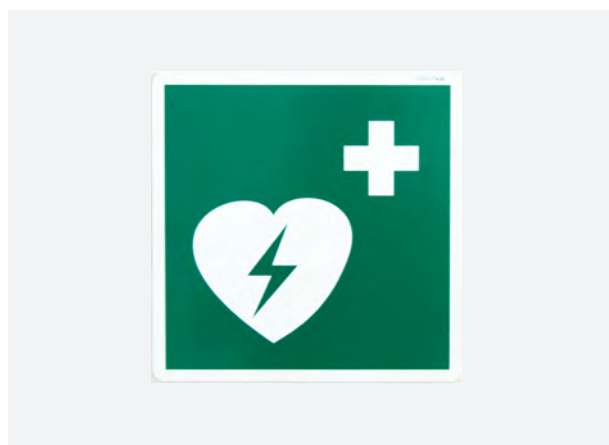

SAFETY SIGN OF AED - ACC170551

Safety Sign AED.
Made of aluminum.
Dimensions: 25x25 cm.

TORCHES

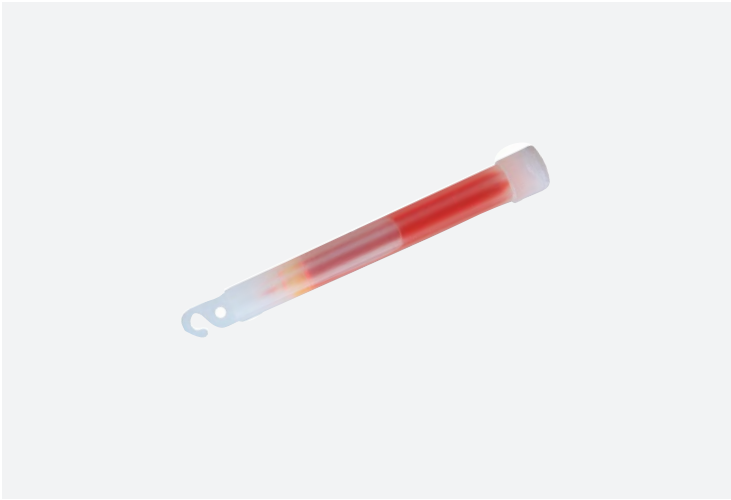

SNAPLIGHT LIGHTSTICK - ACC170044

Bright torch 15 cm long. Ideal emergency light for every occasion. Different usages in the night and in total darkness. Available in brilliant red light. 12 hours duration.



TORCH "SMALL TRAFFIC BATON" 18 LED RED - ACC170113

Fluorescent LED electronic signal torch. Superior to whatever conventional torch concerning high quality and 400 hours autonomy. Waterproof and very resistant to impacts. A magnetic base allowing fixing on metallic surfaces.

Batteries required: 3x1,5 V (not included)

Dimensions: about 28x4 cm

Weight: about 200 gram



TORCH "TRAFFIC LED BATON" 12 LED

Electronic warning torch with blinking light in red or yellow LED. More than any conventional torch for high quality and for the autonomy of 400 hours. Waterproof and very resistant to impacts. Dimensions: about 4x54 cm.

ACC170111 - YELLOW LED

ACC170117 - RED LED



OSCAR BOSCAROL SRL

Via Enzo Ferrari 29, 39100 Bolzano - ITALY

VAT NUMBER: 01458460217

PH. +39 0471932893

FAX +39 0257760140

@ info@boscarol.it

@ info@pec.boscarol.it

WWW. boscarol.it

f [OscarBoscarolSrl](#)

in [oscar-boscarol-srl](#)

WORKING SCHEDULE

Monday - Tuesday

09:00-12:30 and 14:00-18:00

Friday: 09:00-13:00

SELLING CONDITIONS

Selling conditions

1. OVERALL

General selling conditions listed hereupon are to be considered valid for every purchasing made to our company Oscar Boscarol S.r.l. and make a fundamental requirement for every purchase order. Possible variations to such clauses rather than particular conditions must be explicitly agreed in written form.

2. ORDERS

Purchase orders are to be sent to the following: fax 0039 0257760140; email info@boscarol.it or on-line via our web site www.boscarol.it (against registration). All orders received will be processed by our commercial department and must be approved by the Oscar Boscarol S.r.l. through a regular order confirmation. Orders will have to be compliant to the packaging details stated on the price list.

3. TRANSPORT

Transport occurs to the customer's account with the amount stated on the invoice or as assigned port (the customer arrange the shipment with his own forwarder). The customer has to check the packaging status at the moment of delivery: should any anomaly be noticed the goods should be accepted with note to the forwarder itself and to the Oscar Boscarol S.r.l.

4. DELIVERY

Delivery terms are stated on our order confirmations and might be subject to variations which are communicated to the customer. Order processing is strictly subject to its written confirmation from the customer.

5. PRICES

Price listed hereupon are to be considered in Euro, VAT and tax excluded. OSCAR BOSCAROL S.r.l. owns the right to perform changes on the prices in case raw materials and/or components price increase occurred.

6. RETURN POLICY AND COMPLAIN

Possible complains must be sent in written form to Oscar Boscarol S.r.l. within 8 days from goods delivery. The return of the goods must be previously authorized by Oscar Boscarol S.r.l.

The authorized returning goods must be addressed to the Oscar Boscarol S.r.l. warehouse at the conditions stated on the return authorization and needs to be shipped together with a transport document stating the reasons of the return together with purchasing references. Return goods will be replaced and/or credited after check-in inspection.

7. SIZE AND DRAWINGS

Indications related to dimensions and drawings on each publication issued by Oscar Boscarol srl have a mere informative purpose and can be modified without notice.

8. PRODUCTS STOCK AVAILABILITY:

All the products marked with (*) are to be intended as AVAILABLE UNTIL STOCK AVAILABILITY.

9. PAYMENT

Payment must be performed in due course and for the amount agreed. Failing to do so might result in the temporarily suspension of supplies and start bank interests appliance.

10. Refer to Bolzano forum for every controversy that might occur.

Please be advised that all the information featured in this catalogue are merely indicative and subject to continuous changes. For any updated information please take a look at our website www.boscarol.it or contact our commercial department.

Oscar Boscarol company reserves the right to update what stated in the technical and advertising files without any prior notice. Please contact our offices via telephone +39 0471 932893 or via email to double check such information.



OSCAR BOSCAROL SRL
Via Enzo Ferrari 29, 39100 Bolzano ITALY - Ph. +39 0471932893 - Fax +39 0257760140
info@boscarol.it - www.boscarol.it

2021/05-EN